

Compat[®] DualPort

Single Lumen Dual Port Nasogastric Feeding and Drainage Tube

Compat[®] DualPort is a single lumen, dual port nasogastric tube for: 1) the administration of enteral nutrition, fluids and medication via the ENFit[®] enteral-specific feeding port and 2) gastric drainage / intermittent aspiration of stomach contents via the funnel type drainage port.

Features at a glance

Compat[®] DualPort features a funnel type as well as an ENFit[®] connector to help **simplify gastric drainage and administration of nutrition, fluids and medications in the ICU** by the use of one tube for both operations.

- Allows both gastric feeding and drainage/aspiration via a single tube
- Avoids the need for a second tube placement when the patient transitions from requiring drainage/ aspiration to enteral feeding.
- Designed to resist collapsing when checking gastric residuals


Clamps help prevent fluid leakage and **ease** gastric fluid aspiration.

Compat[®] DualPort is made of soft radio-opaque polyurethane with a rounded distal end to help promote **patient comfort**.

A pre-lubricated and open “flow-through” stylet is included to **aid in placement**:

- Allows injecting water and radio-opaque fluids without stylet removal
- May help lower the risk of tube displacement on stylet removal

Compat[®] DualPort is **compatible** with all ENFit[®] administration sets and syringes as well as most drainage/ suction devices*.

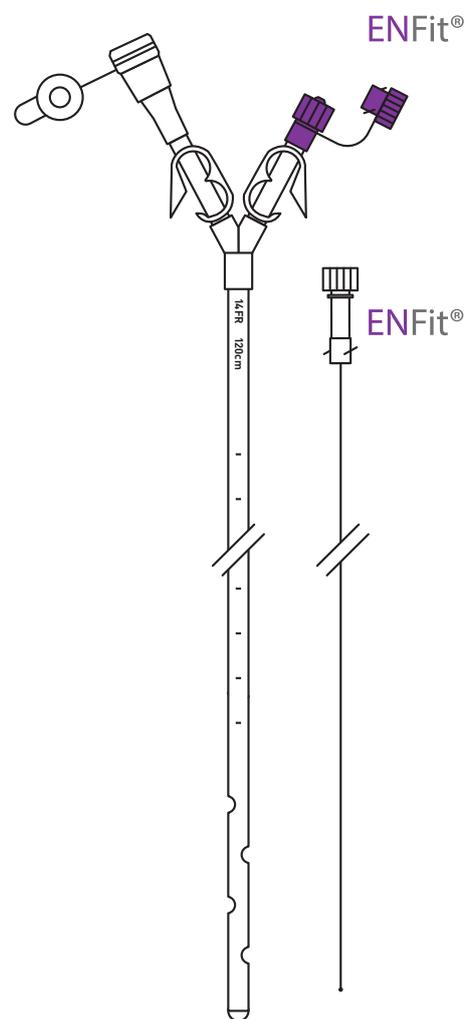


Ordering Information

Product Name	Product Code	Qty/Box
Compat [®] DualPort 14FR 120 CM S	12391311	20

Patient Care. Our Priority.

Single Lumen Dual Port Nasogastric Feeding and Drainage Tube

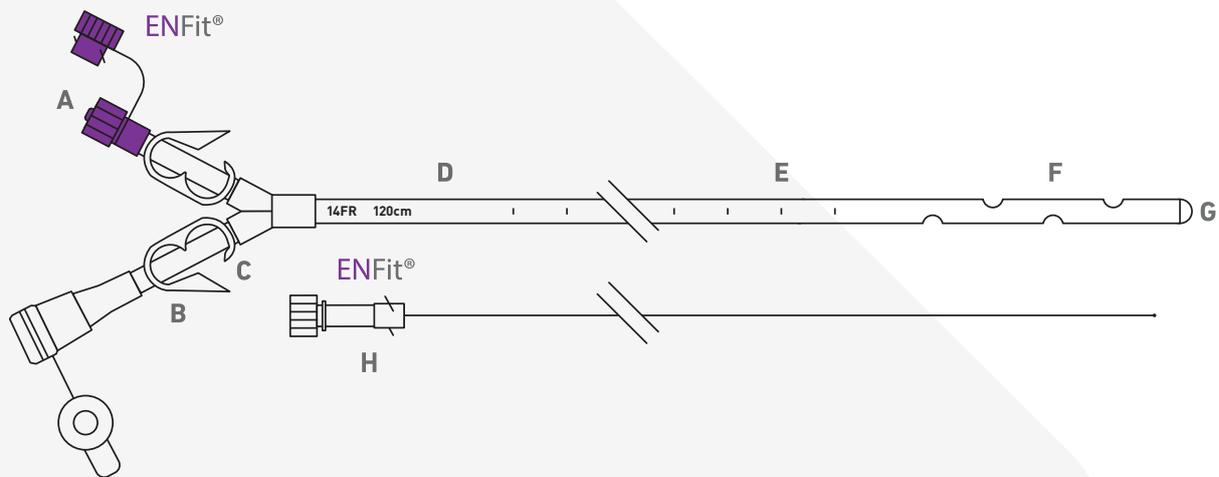

Components

- A ENFit® twist-lockable feeding connector with cap
- B Funnel connector with cap designed to be compatible with most drainage/suction devices*
- C Closing clamps
- D Product identification, printed French (FR) size and length on the shaft
- E Radio-opaque polyurethane tube shaft with placement depth markers every 5 cm
- F 4 lateral alternate openings
- G Rounded distal end
- H Pre-lubricated and open "flow-through" stylet

Specifications

- Sterile:** Sterilized with Ethylene Oxide (ETO).
- Packing:** Individually packed in cases of 20 pieces (29.5 x 25 x 17 cm).
- Storage:** Store in the original packaging at a temperature between +5 and +40°C and in a dry environment. Avoid direct sunlight and heat radiation.
- Shelf life:** 12 to 36 months from production date. Expiry date on the package.
- Duration of use:** Up to 29 days. Single Use.
- Tube length:** 120 cm.
- Tube diameter:** 14FR.
- ISO Standard:** ENFit® connector made according to ISO 80369-3.
- CE Mark:** Complies with the European Medical Device Directive 93/42 EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
- Legal manufacturer:** Cedic Srl, Via Liberazione 63/9, IT-20068 Peschiera Borromeo (MI).
- Ecological footprint:** The 'Grüne Punkt' symbol signals that we contribute to the financing of recycling.
- Materials:** Tubing: Polyurethane. ENFit® connector: Copolyester. Funnel connector: Polyurethane. DEHP, BPA and Latex-free**.
- Country of origin:** Made in Italy.

*Some incompatibilities may exist.
**BPA, Latex and Phthalates including DEHP are not intentionally added during the manufacturing process.
To obtain further information please contact your Compat® representative or go to www.compat.com
*Unless stated otherwise, all trademarks are used with permission.