

# Compat ella<sup>®</sup>

## Compat Ella<sup>®</sup> Dual Spike Set

The Compat Ella<sup>®</sup> Dual Spike Set is intended for the simultaneous administration of enteral nutrition and/or hydration via two separately spiked containers and is suitable for use with the Compat Ella<sup>®</sup> feeding pump.

### Features at a glance

**Dual Spike:** For simultaneous feeding and/or hydration via two separately spiked containers.


**Built-in Safety:** All Compat Ella<sup>®</sup> administration sets feature patented fail-safe **Auto-Anti-Free-Flow (AAFF)** technology. **AAFF** is incorporated into the cassette and is always automatically activated to ensure patient safety.

**Safer Connectors:** All Compat Ella<sup>®</sup> administration sets include ENFit<sup>®</sup> connectors made according to **ISO 80369-3**. **ENFit<sup>®</sup>** and **ENPlus connectors** are designed to be **IV Incompatible** to promote **patient safety** by helping to avoid accidental misconnections.

**Food Safe Materials:** All Compat Ella<sup>®</sup> administration sets are **DEHP, BPA and Latex-Free\***.

**Convenient:**

- Nutrition containers: ENPlus
- Enteral feeding tubes: ENFit<sup>®</sup>
- Medication dispensers: ENFit<sup>®</sup>

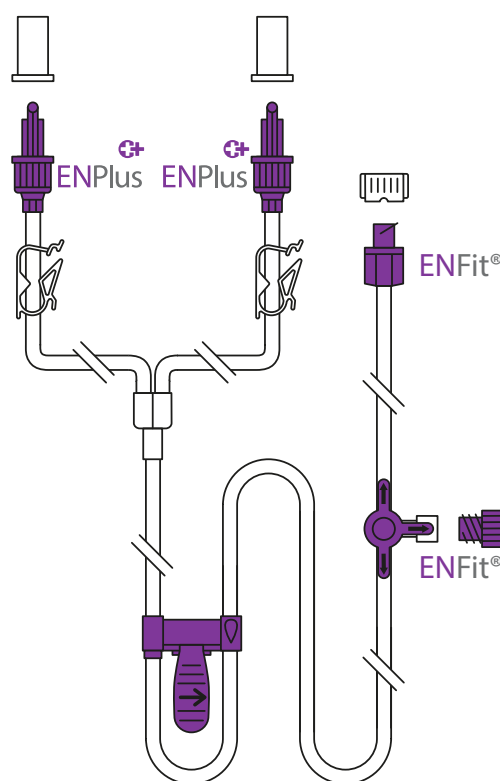


### Ordering Information

Product Name	Product Code	Qty/Box
Compat Ella <sup>®</sup> Dual Spike Set	12297288	30

Have you met Compat Ella<sup>®</sup>?

# Compat Ella<sup>®</sup> Dual Spike Set

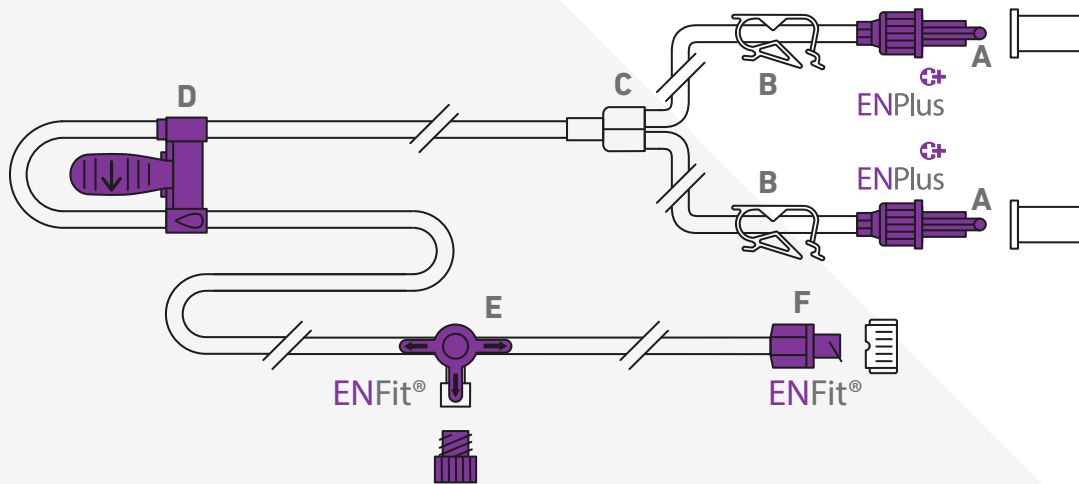

## Components

- A ENPlus 'enteral only' proximal end (nutrition connection) with revolving lock and cap
- B Pinch clamp to manually shut off the flow of fluids
- C Y junction
- D Compat Ella<sup>®</sup> pump cassette with integrated and patented Auto-Anti-Free-Flow (AAFF) safety feature
- E ENFit<sup>®</sup> 3-way medication port with cap
- F ENFit<sup>®</sup> distal end with cap (patient side tube connection)

## Specifications

- Sterile:** Sterilized with Ethylene Oxide.
- Packing:** Individually packed in cases of 30 pieces (29.5 x 26 x 19 cm).
- Storage:** Store at room temperature. Keep away from sunlight and humidity.
- Shelf life:** 3 years from production date. Expiry date on the package.
- Duration of use:** 24 hours. Single Use.
- Tube length:** Distance from nutrition connection to patient side tube connection: 230 cm.  
Distance from nutrition connection to Y junction: 30 cm.  
Distance from nutrition connector to pump: 50 cm.
- Tube diameters:** Nominal inner diameter: 3.0 mm. Nominal outer diameter: 4.15 mm.
- CE mark:** Complies with the European Medical Device Directive 93/42 EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
- Legal manufacturer:** Cedic Srl, Via Liberazione 63/9, IT-20068 Peschiera Borromeo (MI).
- Ecological footprint:** The 'Grüne Punkt' symbol signals that we contribute to the financing of the recycling.
- Materials:** Tubing: PVC and silicone. All materials DEHP, BPA and Latex-free\*. Additional details on request.
- Country of origin:** Made in Italy.

\*Some incompatibilities may exist.  
\*\*BPA, Latex and Phthalates including DEHP are not intentionally added during the manufacturing process.

To obtain further information please contact your Compat<sup>®</sup> representative or go to [www.compat.com](http://www.compat.com)  
\*Unless stated otherwise, all trademarks are used with permission.

AXI00404 REV.0 14 May 2019