

# Compat<sup>®</sup> Extension Set

## Low-Profile Gastrostomy Tube Extension Sets

The Compat<sup>®</sup> Low-Profile Gastrostomy Tube Extension Sets are intended for the administration of enteral nutrition, fluids and medication. Indicated for patients who require nutritional support and have a low-profile gastrostomy tube (button).

### Features at a glance

#### Compatibility:


- Compat<sup>®</sup> Extension Sets are compatible with most low-profile gastrostomy tubes including the Mic-Key<sup>®</sup> Buttons, the Freka<sup>®</sup> Belly Buttons and the AMT MiniONE<sup>®</sup> Buttons.
- Only compatible with administration sets and syringes with ENFit<sup>®</sup> connectors.

#### Multiple sizes and features:

- Right-angle connector for continuous feeding and straight connector for bolus feeding.
- 5, 30 and 60 cm to address multiple patient and clinician needs.
- Y port to facilitate medication administration.
- Screw cap and closing clamp to help prevent fluid leakage and ease handling.

**Patient safety:** ENFit<sup>®</sup> connectors are made according to ISO 80369-3. They are designed to be IV INcompatible to promote patient safety by helping to avoid accidental misconnections.

**Safe Materials:** All Compat<sup>®</sup> Extension Sets are DEHP, BPA and Latex-Free\*.

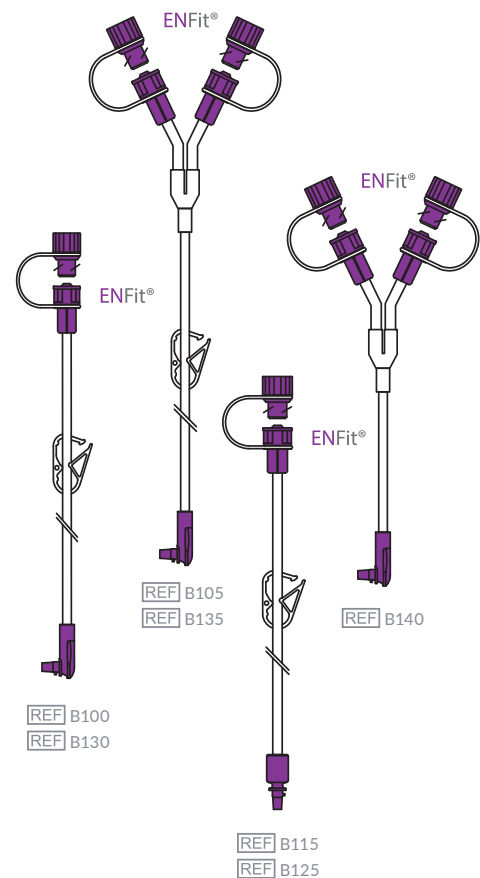


### Ordering Information

	ENFit <sup>®</sup> Connector	Length (cm)	Product Code	Qty/box
Compat <sup>®</sup> ENFit <sup>®</sup> Extension Set with right-angle connector	Dual	5	B140	10
Compat <sup>®</sup> ENFit <sup>®</sup> Extension Set with right-angle connector	Dual	30	B105	
Compat <sup>®</sup> ENFit <sup>®</sup> Extension Set with right-angle connector	Dual	60	B135	
Compat <sup>®</sup> ENFit <sup>®</sup> Extension Set with right-angle connector	Single	30	B100	
Compat <sup>®</sup> ENFit <sup>®</sup> Extension Set with right-angle connector	Single	60	B130	
Compat <sup>®</sup> ENFit <sup>®</sup> Extension Set with straight connector	Single	30	B115	
Compat <sup>®</sup> ENFit <sup>®</sup> Extension Set with straight connector	Single	60	B125	

Patient Care. Our Priority.

# Low-Profile Gastrostomy Tube Extension Sets

ENFit®

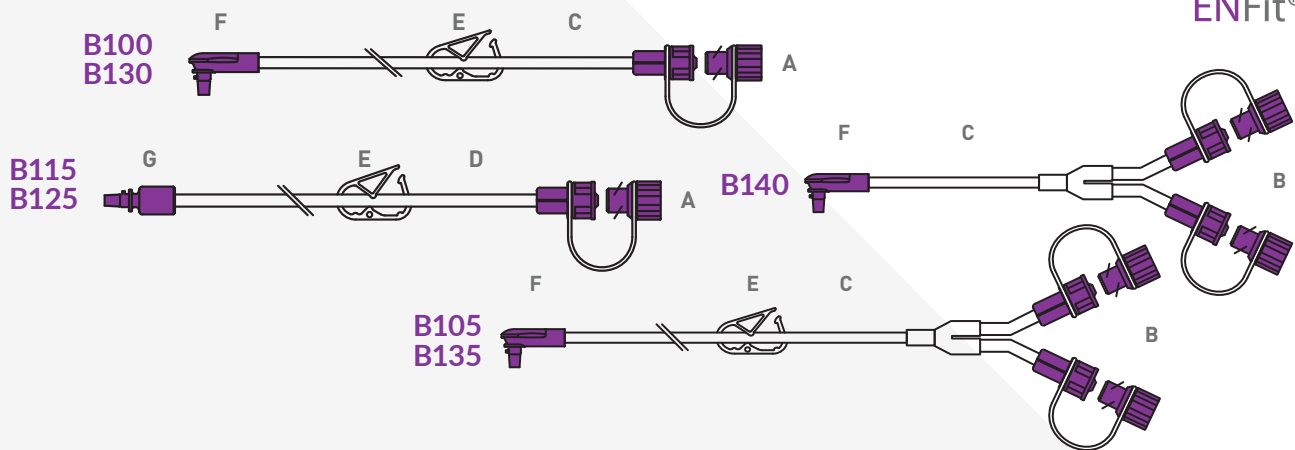

## Components

- A** ENFit® twist-lockable feeding connector with cap and cap tether
- B** Y ENFit® twist-lockable feeding port with cap and tether
- C** Standard PVC tubing 4.15 mm
- D** Large bore PVC tubing 5.4 mm
- E** Closing clamp
- F** Right-angle connector with alignment marker (patient side tube connection)
- G** Straight connector with alignment marker (patient side tube connection)

## Specifications

- Sterile:** Sterilized with Ethylene Oxide.
- Packing:** Individually packed in a box of:  
10 pieces per case (19 x 19 x 10 cm).  
24 cases in a shipping box (60 x 40 x 42 cm).
- Storage:** Store dry at room temperature and away from sunlight.
- Shelf life:** 3 years from production date. Expiry date on the package.
- Duration of use:** Up to 7 days and 5 usages per day.
- Tube length:** 5, 30 or 60 cm.
- Tube diameter:** 4.15 mm (REF with right-angle connector) or 5.4 mm (REF with straight connector).
- ISO Standard:** ENFit® connector made according to ISO 80369-3.
- CE Mark:** Complies with the European Medical Devices Directive 93/42/EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
- Legal manufacturer:** Cedic Srl, Via Liberazione 63/9, IT-20068 Peschiera Borromeo (MI).
- Sustainability:** Environmental management system complies with ISO 14001.
- Materials:** Tubing: PVC. ENFit® and button connectors: Copolyester. All materials DEHP, BPA and Latex-free\*. Additional details on request.
- Country of origin:** Made in Italy.

\*Does not contain phthalates, including DEHP. Not made with natural rubber latex. Not formulated with BPA.  
ENFit® and Compat® trademarks are used with permission. Mic-Key® is a registered trademark of Halyard Health, Inc. or its affiliates.  
Freka® is a registered trademark of Fresenius Kabi AG. MiniONE® is a registered trademark of Applied Medical Technology, Inc.

To obtain further information please contact your Compat® representative or go to [www.compat.com](http://www.compat.com)

AXI00013 REV.1 18 JUNE 2020

# Compat®