

# Compat<sup>®</sup> Flexibaggle<sup>®</sup>

## Foldable Empty Enteral Reservoir 1.5L

The Compat<sup>®</sup> Flexibaggle<sup>®</sup> 1.5L foldable empty reservoir with ENPlus port is intended for the administration of enteral nutrition and / or fluids to patients.

### Features at a glance


The Compat<sup>®</sup> Flexibaggle<sup>®</sup> 1.5L foldable empty reservoir with ENPlus port has been designed to help improve patient care.

**Ease of use:** The built-in hanger and large top-fill opening allow the bag to be easily filled with nutrition or water and its fold-ability saves on storage space enabling it to easily be at hand when needed.

**Ambulatory use:** The bag resists a fall of up to 1.5m and the easy-to-close cap design requires only a quarter turn to ensure a secure seal, allowing the product to also be used on-the-go with a low risk of leakage.

**Patient safety:** Compat<sup>®</sup> Flexibaggle<sup>®</sup> is delivered sterile and includes an ENPlus connector made according to ISO 18250-3. The ENPlus connector is designed to be IV INcompatible to promote patient safety by helping to avoid accidental misconnections. Compat<sup>®</sup> Flexibaggle<sup>®</sup> is made of Ethylene-vinyl Acetate (EVA), a soft and food friendly material, and is also DEHP, BPA and Latex-Free\*.

**Eco-friendly:** Compat<sup>®</sup> Flexibaggle<sup>®</sup> is foldable, reducing the environmental impact of shipping.

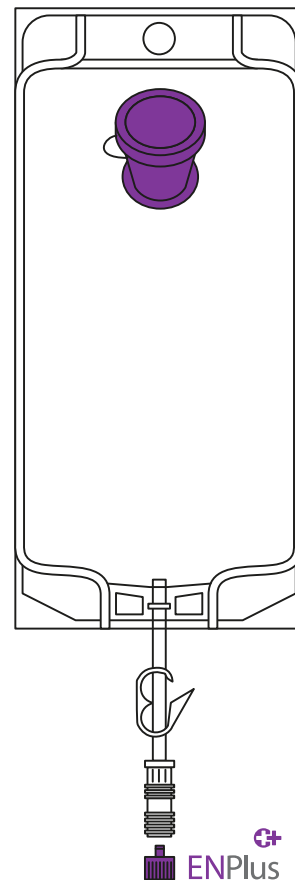


### Ordering Information

Product Name	Product Code	Qty/Box
Compat <sup>®</sup> Flexibaggle <sup>®</sup> 1.5L	12306162	30

*Patient Care. Our Priority.*

# Foldable Empty Enteral Reservoir 1.5L

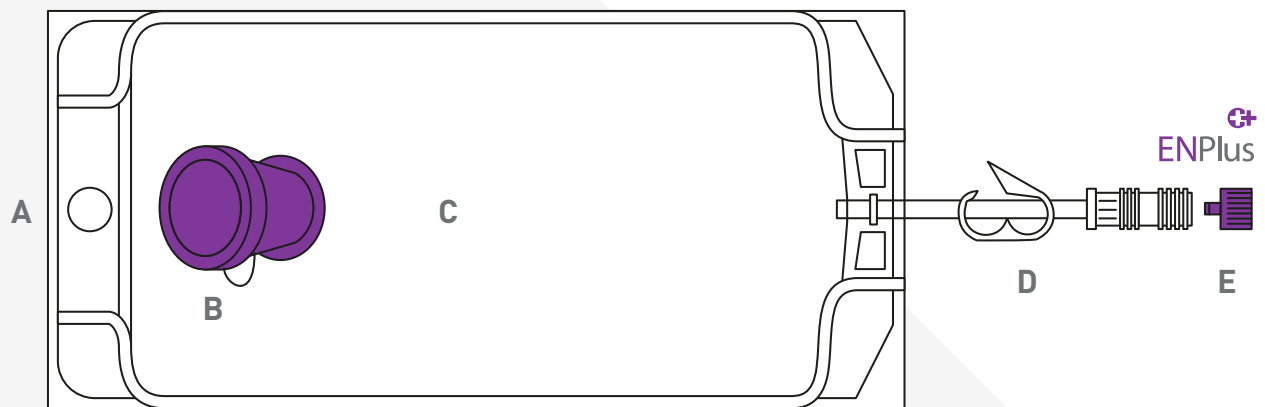

## Components

- A Built-in hanger
- B Large top-fill opening with secure screw cap
- C 1.5L container with graduations (at 50 mL and then every 100 mL) and tags for patient name and other important information (formula, volume, date and rate)
- D Pinch clamp
- E ENPlus 'enteral only' feeding port with cap

## Specifications

- Wide-mouth port:** Nominal inner diameter:  $\pm 38$  mm.
- Bag size:** Max. nominal bag size:  $\pm 336 \times 188$  mm.
- Sterile:** Sterilized with Ethylene Oxide.
- Packing:** Individually packed in cases of 30 pieces (58.5 x 38.5 x 16 cm).
- Storage:** Store dry at room temperature and away from sunlight.
- Shelf life:** Up to 3 years from production date. Expiry date on the package.
- Duration of use:** 24 hours. Single Use.
- ISO Standard:** ENPlus connector made according to ISO 18250-3..
- CE Mark:** Complies with the European Medical Devices Directive 93/42/EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
- Legal manufacturer:** Cedic Srl, Via Liberazione 63/9, IT-20068 Peschiera Borromeo (MI).
- Sustainability:** Environmental management system complies with ISO 14001:2015.
- Materials:** Bag: EVA. Tubing: PVC. All materials DEHP, BPA and Latex-free\*. Additional details on request.
- Country of origin:** Made in Italy.

\*BPA, Latex and Phthalates including DEHP are not intentionally added during the manufacturing process.

To obtain further information please contact your Compat® representative or go to [www.compat.com](http://www.compat.com)  
\*Unless stated otherwise, all trademarks are used with permission.

AXI00032 REV.0 27 NOVEMBER 2019

# Compat®