

Compat[®] Adapter

Compat[®] Adapter Stepped Set > ENFit[®] Tube

The Compat[®] Adapter Stepped Set > ENFit[®] Tube is intended to connect an ENFit[®] feeding tube with a stepped administration set, a large catheter tip syringe or a gastric drainage/suction device.

Features at a Glance

The Compat[®] Adapter Stepped Set > ENFit[®] Tube is placed between the feeding tube and the administration set, the syringe or the gastric drainage/suction device.


Ease of Use: The Compat[®] Adapter Stepped Set > ENFit[®] Tube is designed to **facilitate the transition to ENFit[®]** by allowing compatibility between ENFit[®] feeding tubes and most administration sets equipped with a stepped distal end. It is also compatible with most drainage or aspiration systems **to ease gastric residual volume management**.

Safer Connectors: The Compat[®] Adapter Stepped Set > ENFit[®] Tube includes an ENFit[®] connector made according to ISO 80369-3. ENFit[®] is designed to be IV incompatible in order to help **promote patient safety** by helping to avoid accidental misconnections.

Safe Materials: The Compat[®] Adapter Stepped Set > ENFit[®] Tube is **DEHP, BPA and Latex-Free***.

Compatibility:

- Enteral feeding tubes: ENFit[®]
- Administration Sets: Stepped, ENLock
- Syringes: Large Catheter Tip, ENLock
- Drainage/suction devices: Funnel**

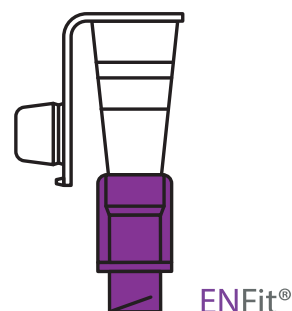


Ordering Information

Product Name	Product Code	Qty/Box
Compat [®] Adapter Stepped Set > ENFit [®] Tube	12319560	30 pcs

Patient Care. Our Priority.

Compat® Adapter Stepped Set > ENFit® Tube

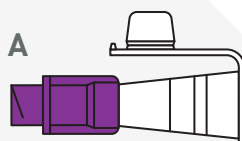

Components

A Adapter Stepped Set with cap

Specifications

Sterilization:	Non-sterile product.
Packing:	Individually packed in cases of 30 pieces (19 × 19 × 10 cm).
Storage:	Store at room temperature. Keep away from sunlight and humidity.
Shelf life:	3 years from production date. Expiry date on the package.
Duration of use:	24 hours. Single Use.
ISO Standard:	ENFit® connector made according to ISO 80369-3.
CE Mark:	Complies with the European Medical Device Regulation 2017/745.
Classification:	Class I.
Legal manufacturer:	Cedic Srl, Via Liberazione 63/9, IT-20068 Peschiera Borromeo (MI).
Sustainability:	Environmental management system complies with ISO 14001:2015.
Materials:	All materials DEHP, BPA and Latex-free*. Additional details on request.
Country of origin:	Made in Italy.

*DEHP: Not intentionally added during the manufacturing process. Not formulated with BPA or natural rubber latex.

**Some incompatibilities may exist.

To obtain further information please contact your Compat® representative or go to www.compat.com

*Unless stated otherwise, all trademarks are used with permission.

This document does not replace the original Instructions for Use which must be read before use.

AXI01534 Rev0 24 NOVEMBER 2022

Compat®