

Compat ella[®]

Compat Ella[®] Universal Set

The Compat Ella[®] Universal Set is intended for administration of enteral nutrition, fluids or medication and is suitable for use with the Compat Ella[®] feeding pump.

Features at a glance

The Compat Ella[®] Universal Set features an ENPlus nutrition connector, a universal screw cap with air filter, a 3-way medication port and an ENFit[®] distal end (patient side tube connection). A separate bottle hanger is also included.


Built-in Safety: All Compat Ella[®] administration sets feature patented fail-safe Auto-Anti-Free-Flow (AAFF) technology. AAFF is incorporated into the cassette and is always automatically activated to ensure patient safety.

Safer Connectors: All Compat Ella[®] administration sets include ENFit[®] connectors made according to ISO 80369-3. The ENFit[®] and ENPlus connectors are designed to be IV Incompatible to promote patient safety by helping to avoid accidental misconnections.

Food Safe Materials: All Compat Ella[®] administration sets are DEHP, BPA and Latex-Free*.

Compatibility:

- Nutrition containers: ENPlus, glass and plastic wide-neck containers, crown cork glass bottles
- Enteral feeding tubes: ENFit[®]
- Medication dispensers: ENFit[®]

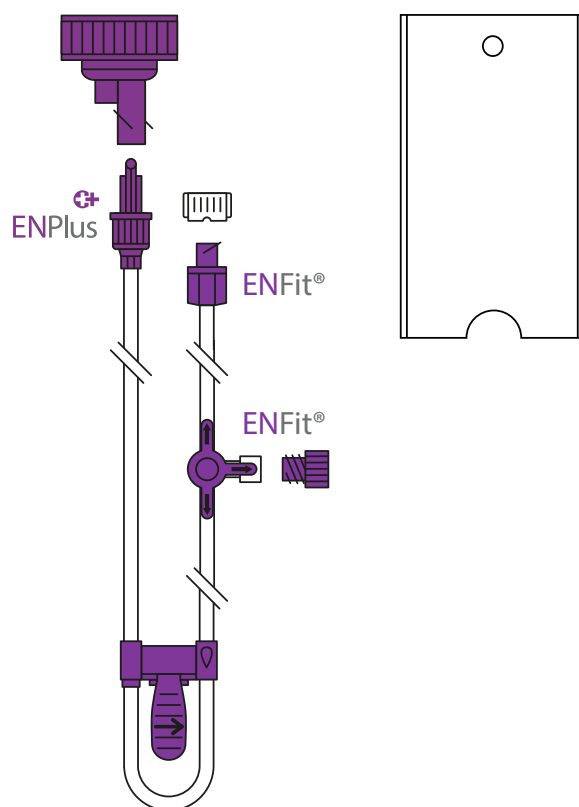


Ordering Information

Product Name	Product Code	Qty/Box
Compat Ella [®] Universal Set	12298206	30

Have you met Compat Ella[®]?

Compat Ella® Universal Set

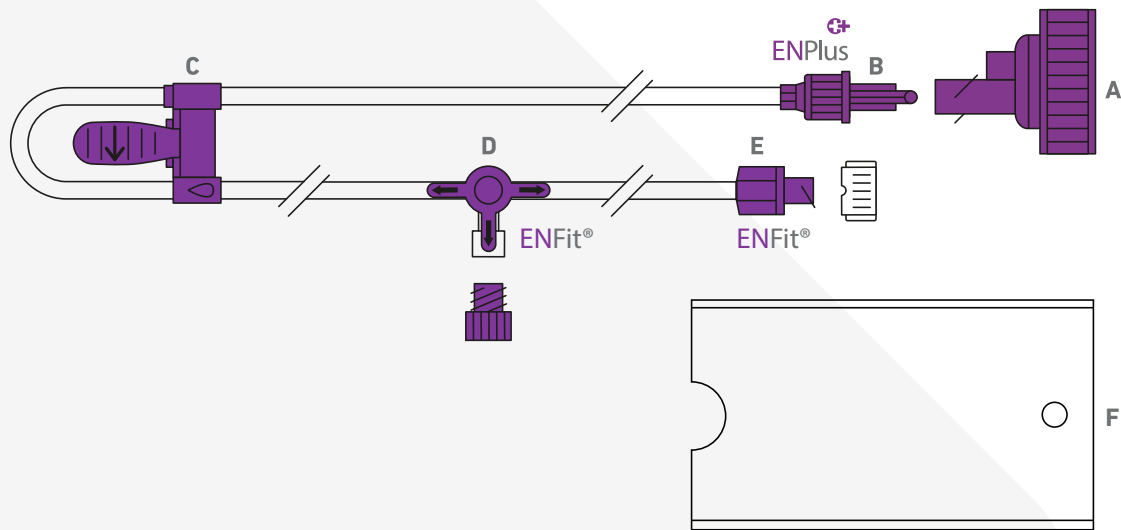

Components

- A Universal screw cap with air filter and ENPlus port for use with glass and plastic wide-neck containers (DIN 168-1) or crown cork glass bottles (EN 14634 and EN 14635)
- B ENPlus 'enteral only' proximal end (nutrition connection) with revolving lock
- C Compat Ella® pump cassette with integrated and patented Auto-Anti-Free-Flow (AAFF) safety feature
- D ENFit® 3-way medication port with cap
- E ENFit® distal end with cap (patient side tube connection)
- F 500 mL bottle holder LDPE (Low Density Polyethylene)

Specifications

Sterile:	Sterilized with Ethylene Oxide.
Packing:	Individually packed in cases of 30 pieces (39.5 x 29.5 x 14 cm).
Storage:	Store at room temperature. Keep away from sunlight and humidity.
Shelf life:	3 years from production date. Expiry date on the package.
Duration of use:	24 hours. Single Use.
Tube Length:	Distance from nutrition connection to patient side tube connection: 230 cm. Distance from nutrition connector to pump: 50 cm.
Tube diameters:	Nominal inner diameter: 3.0 mm. Nominal outer diameter: 4.15 mm.
CE mark:	Complies with the European Medical Device Directive 93/42 EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
Legal manufacturer:	Cedic Srl, Via Liberazione 63/9, IT-20068 Peschiera Borromeo (MI).
Ecological footprint:	The 'Grüne Punkt' symbol signals that we contribute to the financing of the recycling.
Materials:	Tubing: PVC and silicone. All materials DEHP, BPA and Latex-free*. Additional details on request.
Country of origin:	Made in Italy.

*BPA, Latex and Phthalates including DEHP are not intentionally added during the manufacturing process.

To obtain further information please contact your Compat® representative or go to www.compat.com

®Unless stated otherwise, all trademarks are used with permission.

AXI00406 REV.0 15 JUNE 2020

Compat®