

# Compat<sup>®</sup> Buddy

## Full Range of Low-profile Gastrostomy Feeding Tube Extension Sets

The Compat<sup>®</sup> Buddy Low-profile Gastrostomy Feeding Tube Extension Sets are intended for the administration of enteral nutrition, fluids and/or medication. Indicated for patients who require nutritional support and have a low-profile gastrostomy feeding tube (button). Also suitable for passive gastric decompression.

### Features at a Glance

#### Multiple Sizes and Features:


- Right-angle connector for continuous feeding and straight connector for bolus feeding.
- 5, 30 and 60 cm to address multiple patient and clinician needs.
- Dual port to facilitate medication administration (optional).
- Screw cap and closing clamp (optional) to help prevent fluid leakage and ease handling.

#### Compatibility:

- Compat<sup>®</sup> Buddy Extension Sets are compatible with Compat<sup>®</sup> Buddy button and with most other buttons including Mic-Key<sup>®</sup>, Aquarius<sup>™</sup>, MiniONE<sup>®</sup> and Freka<sup>®</sup> Belly.
- Only compatible with administration sets and syringes with ENFit<sup>®</sup> connector to promote patient safety.

**Patient Safety:** ENFit<sup>®</sup> connectors are made according to ISO 80369-3. They are designed to be IV incompatible to promote patient safety by helping to avoid accidental misconnections.

**Safe Materials:** DEHP, BPA and Latex-Free\*.

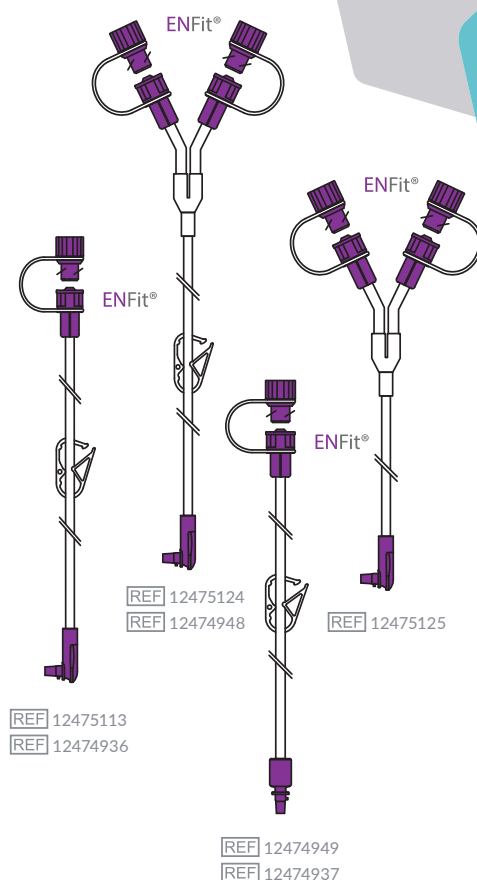


### Ordering Information

	Connector	ENFit <sup>®</sup> Port	Length (cm)	Product Code	Qty/box
Compat <sup>®</sup> Buddy Low-profile Gastrostomy Feeding Tube Extension Set	Right-angle	Dual	5	12475125	10
	Right-angle	Dual	30	12475124	
	Right-angle	Dual	60	12474948	
	Right-angle	Single	30	12475113	
	Right-angle	Single	60	12474936	
	Straight	Single	30	12474949	
	Straight	Single	60	12474937	

*Patient Care. Our Priority.*

# Full Range of Low-profile Gastrostomy Feeding Tube Extension Sets

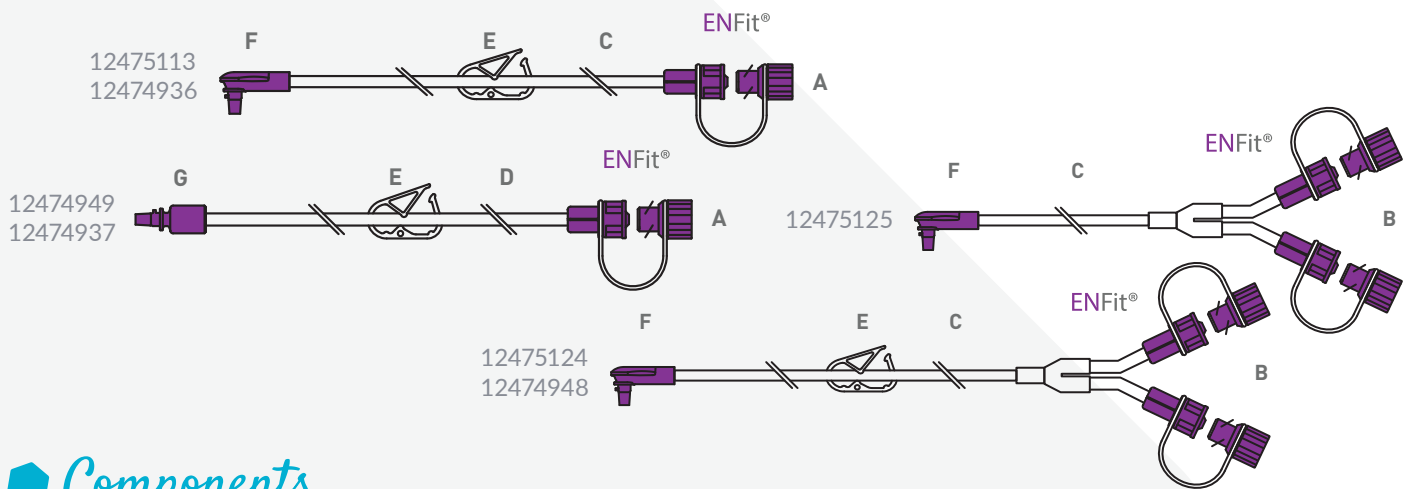

## Components

- A** ENFit® twist-lockable feeding port with cap and tether
- B** Dual ENFit® twist-lockable feeding port with cap and tether
- C** Standard PVC tubing outer diameter 4.15 mm
- D** Large bore PVC tubing outer diameter 5.4 mm
- E** Closing clamp
- F** Right-angle connector with alignment marker (patient side tube connection)
- G** Straight connector with alignment marker (patient side tube connection)

## Specifications

- Sterile:** Sterilized with Ethylene Oxide (EO).
- Packing:** Individually packed in a box of:  
10 pieces per case (19 x 19 x 10 cm).  
24 cases in a shipping box (60 x 40 x 42 cm).
- Storage:** Store dry at room temperature and away from sunlight.
- Shelf life:** 3 years from production date. Expiry date on the package.
- Duration of use:** Up to 7 days and 5 usages per day.
- Tube length:** 5, 30 or 60 cm.
- Outer tube diameter:** 4.15 mm (REF with right-angle connector) or 5.4 mm (REF with straight connector).
- ISO Standard:** ENFit® connector made according to ISO 80369-3.
- CE mark:** Complies with the European Medical Devices Directive 93/42/EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
- Classification:** Class IIa.
- Legal manufacturer:** Cedic Srl, Via Liberazione 63/9, IT-20068 Peschiera Borromeo (MI).
- Sustainability:** Environmental management system complies with ISO 14001:2015.
- Materials:** Tubing: PVC. ENFit® and button connectors: Copolyester. All materials DEHP, BPA and Latex-free\*. Additional details on request.
- Country of origin:** Made in Italy.

\*DEHP and BPA are not intentionally added during the manufacturing process. Not formulated with natural rubber latex. ENFit® and Compat® trademarks are used with permission. Mic-Key® is a registered trademark of Halyard Health, Inc. or its affiliates. Freka® is a registered trademark of Fresenius Kabi AG. MiniONE® is a registered trademark of Applied Medical Technology, Aquarius™ is a registered trademark of Degania Medical, Inc.

To obtain further information please contact your Compat® representative or go to [www.compat.com](http://www.compat.com)

AXI00469 REV.3 16 MARCH 2021