

# Compat<sup>®</sup> Buddy

## Stoma Measuring Device

Compat<sup>®</sup> Buddy stoma measuring device is designed to determine the depth of a well-established gastrostomy stoma tract in order to assist in proper length selection of the Compat<sup>®</sup> Buddy gastrostomy button.


### Features at a Glance

Compat<sup>®</sup> Buddy stoma measuring device features a tubular shaft with graduated scale, a movable disk and a retention balloon to measure the length of the stoma.

**Patient Comfort:** Compat<sup>®</sup> Buddy stoma measuring device allows to select the proper length of the Compat<sup>®</sup> Buddy gastrostomy button to help prevent gastric leakage and skin irritation.

**Safe Materials:** Compat<sup>®</sup> Buddy stoma measuring device is sterile, DEHP, BPA and Latex-free\*.

**A Range for Most Needs:** Compat<sup>®</sup> Buddy stoma measuring device belongs to a full range of gastrostomy buttons and accessories designed to meet most clinical needs for long-term enteral feeding.

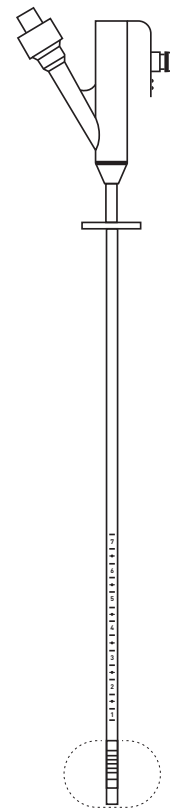

### Ordering Information

Product Name	Product Code	Qty/box
Compat <sup>®</sup> Buddy Stoma Measuring Device	12461149	2

*Patient Care. Our Priority.*

# Stoma Measuring Device



## Components

- A Inflation valve with Luer connector
- B Funnel connector with closing cap
- C Movable disk
- D 12FR tube shaft with graduated scale every centimeter and markers every 0.25 cm
- E 5 ml silicone retention ballon

## Specifications

- Sterile:** Sterilized with Ethylene Oxide (EO).
- Packaging:** Individually packed in a box of 2 pieces.
- Storage:** Store in the original packaging at room temperature and in a dry environment. Avoid direct sunlight and heat radiation.
- Shelf life:** 3 years from production date. Expiry date on the package.
- Duration of use:** For single use.
- Size:** 12FR 5 ml.
- CE Mark:** Complies with the European Medical Devices Directive 93/42/EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
- Classification:** Class I (sterile).
- Legal manufacturer:** Degania Silicone Ltd.
- Materials:** Silicone. All materials are DEHP, BPA and Latex-free\*. Additional details on request.

\*DEHP and BPA are not intentionally added during the manufacturing process. Not formulated with natural rubber latex.

To obtain further information please contact your Compat® representative or go to [www.compat.com](http://www.compat.com)  
\*Unless stated otherwise, all trademarks are used with permission.

AXI00333 REV.1 01 MARCH 2021

# Compat®