

Compat[®] DualPort

Single Lumen Dual Port Nasogastric Feeding and Drainage Tube

Compat[®] DualPort is a single lumen, dual port nasogastric tube for: 1) the administration of enteral nutrition, fluids and medication via the ENFit[®] enteral-specific feeding port and 2) gastric drainage/intermittent aspiration of stomach contents via the funnel type drainage port.

Features at a Glance

Designed for the ICU: Compat[®] DualPort features a funnel type as well as an ENFit[®] connector to help simplify gastric drainage and administration of nutrition, fluids and medication in the ICU by the use of one tube for both operations.

- Allows both gastric feeding and drainage/aspiration via a single tube.
- Prevents the need for a second tube placement when the patient transitions from requiring drainage/aspiration to enteral feeding.
- Designed to resist collapsing when checking gastric residuals.


Clamps help prevent fluid leakage and ease gastric fluid aspiration.

Patient Comfort: Compat[®] DualPort is made of soft radio-opaque polyurethane with a rounded distal end to help promote patient comfort.

Aid in Placement: A pre-lubricated and open “flow-through” stylet is included to aid in placement:

- Allows injecting water and radio-opaque fluids without stylet removal.
- May help lower the risk of tube displacement on stylet removal.

Compat[®] DualPort is compatible with all ENFit[®] administration sets and syringes as well as most drainage/suction devices*.

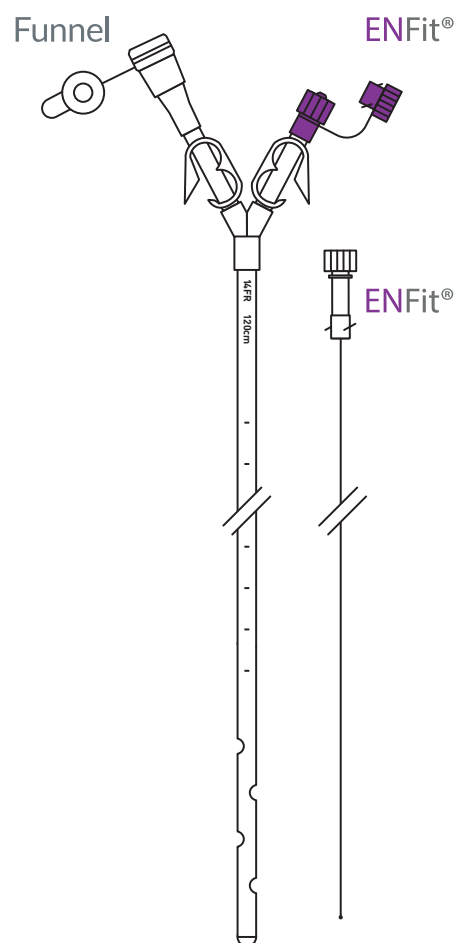


Ordering Information

Product Name	Product Code	Qty/Box
Compat [®] DualPort 14FR 120 CM S	12391311	20

Patient Care. Our Priority.

Single Lumen Dual Port Nasogastric Feeding and Drainage Tube

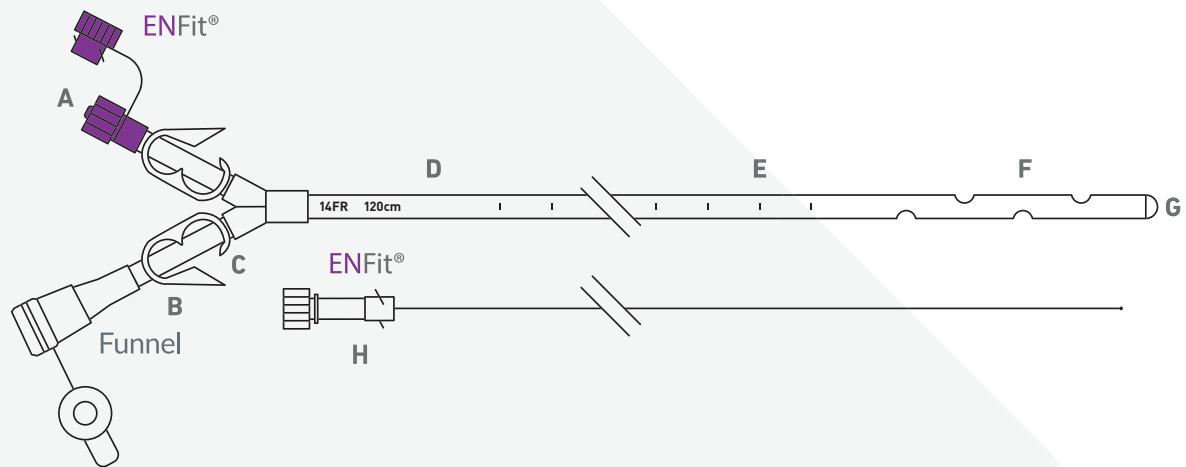

Components

- A ENFit® twist-lockable feeding connector with cap
- B Funnel connector with cap designed to be compatible with most drainage/suction devices*
- C Closing clamps
- D Product identification, printed French (FR) size and length on the shaft
- E Radio-opaque polyurethane tube shaft with placement depth markers every 5 cm
- F 4 lateral alternate openings
- G Rounded distal end
- H Pre-lubricated and open “flow-through” stylet

Specifications

- Sterile:** Sterilized with Ethylene Oxide (EO).
- Packing:** Individually packed in cases of 20 pieces (29.5 x 25 x 17 cm).
- Storage:** Store in the original packaging at room temperature and in a dry environment. Avoid direct sunlight and heat radiation.
- Shelf life:** 3 years from production date. Expiry date on the package.
- Duration of use:** Up to 29 days. Single Use.
- Length:** 120 cm.
- Diameter:** 14FR.
- ISO Standard:** ENFit® connector made according to ISO 80369-3.
- CE Mark:** Complies with the European Medical Devices Directive 93/42/EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
- Classification:** Class IIa.
- Legal manufacturer:** Cedic Srl, Via Liberazione 63/9, IT-20068 Peschiera Borromeo (MI).
- Sustainability:** Environmental management system complies with ISO 14001:2015.
- Materials:** Tubing: Polyurethane. ENFit® connector: Copolyester. Funnel connector: Polyurethane. DEHP, BPA and Latex-free**. Additional details on request.
- Country of origin:** Made in Italy.

*Some incompatibilities may exist.
**DEHP and BPA are not intentionally added during the manufacturing process. Not formulated with natural rubber latex.
To obtain further information please contact your Compat® representative or go to www.compat.com
*Unless stated otherwise, all trademarks are used with permission.