

# Compat ella®

## Compat Ella® Combiset 1.5 L

The Compat Ella® Combiset 1.5 L is intended for administration of enteral nutrition, fluids or medication and is suitable for use with the Compat Ella® feeding pump.

### Features at a Glance

The Compat Ella® Combiset 1.5 L features an integrated 1.5 L Flexibaggle® reservoir, an ENFit® 3-way medication port and an ENFit® distal end (patient side tube connection).

**Reliable Feeding:** The Compat Ella® Combiset 1.5 L features patented fail-safe Auto-Anti-Free-Flow (AAFF) technology. AAFF is incorporated into the cassette and is always automatically activated to ensure patient safety.


**Safer Connectors:** The Compat Ella® Combiset 1.5 L includes ENFit® connectors made according to ISO 80369-3. The ENFit® connector is designed to be IV incompatible to promote patient safety by helping to avoid accidental misconnections.

**Convenience:** The integrated 1.5 L reservoir offers a convenient open system with a built-in hanger and screw cap. The large top-fill opening allows the bag to be refilled with nutrition or water without stopping the feeding program, or removing the reservoir from the hanger.

**Safe Materials:** The 1.5 L reservoir is made of Ethylene-vinyl Acetate (EVA), a soft and food compatible material. The Compat Ella® Combiset 1.5 L is sterile, DEHP, BPA and Latex-Free\*.

**Compatibility:**

- Enteral feeding tubes: ENFit®
- Syringes: ENFit®

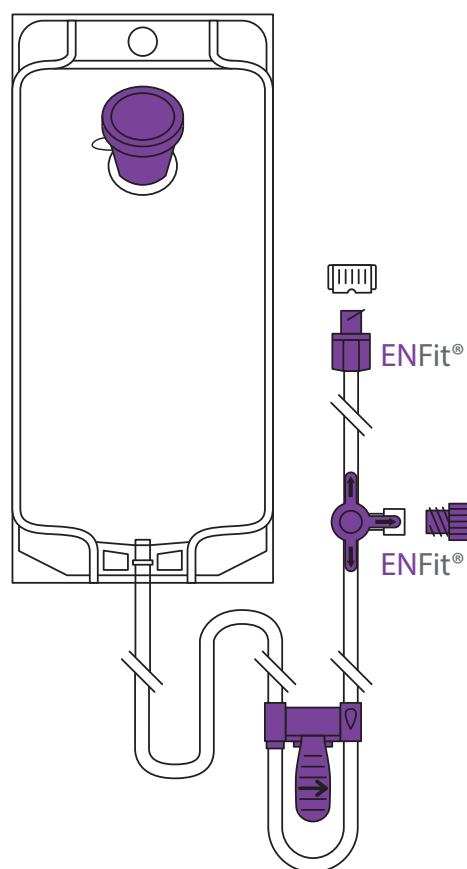


### Ordering Information

Product Name	Product Code	Qty/Box
Compat Ella® Combiset 1.5 L	12298202	30

Have you met Compat Ella®?

# Compat Ella® Combiset 1.5 L

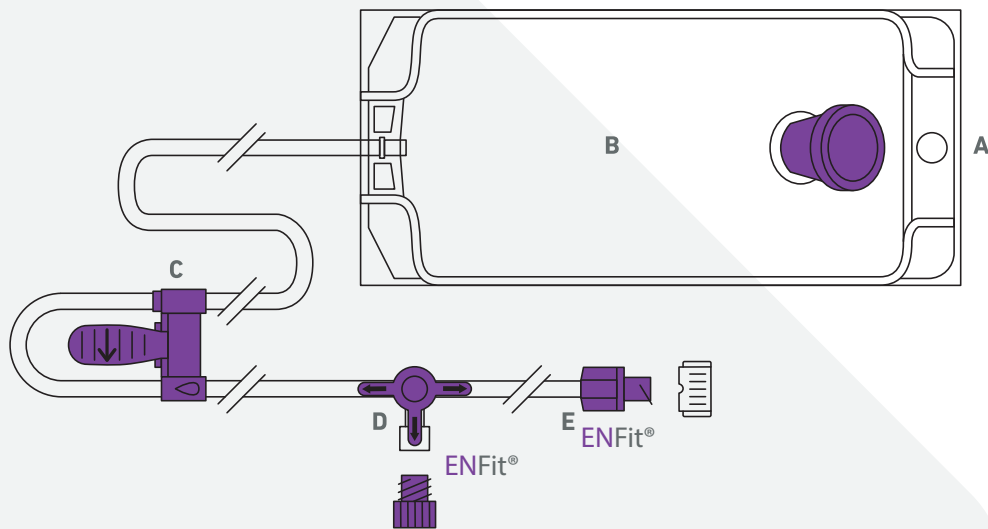

## Components

- A** Large top-fill opening with built-in hanger and screw cap
- B** Flexibaggle® 1.5 L container made of EVA with graduations (at 50 ml and every 100 ml) and tags for patient name and other important information (formula, volume, date and rate)
- C** Compat Ella® pump cassette with integrated and patented Auto-Anti-Free-Flow (AAFF) safety feature
- D** ENFit® compatible 3-way medication port with cap
- E** ENFit® distal end with cap (patient side tube connection)

## Specifications

<b>Sterile:</b>	Sterilized with Ethylene Oxide (EO).
<b>Packing:</b>	Individually packed in cases of 30 pieces (58.5 x 38.5 x 16 cm).
<b>Storage:</b>	Store at room temperature. Keep away from sunlight and humidity.
<b>Shelf life:</b>	3 years from production date. Expiry date on the package.
<b>Duration of use:</b>	24 hours. Single Use.
<b>Length:</b>	Distance from bag to patient side tube connection: 230 cm. Distance from Flexibaggle® to pump: 50 cm.
<b>Diameters:</b>	Nominal inner diameter: 3.0 mm. Nominal outer diameter: 4.15 mm.
<b>ISO Standard:</b>	ENFit® connector made according to ISO 80369-3.
<b>CE Mark:</b>	Complies with the European Medical Devices Directive 93/42/EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
<b>Classification:</b>	Class IIa.
<b>Legal manufacturer:</b>	Cedic Srl, Via Liberazione 63/9, IT-20068 Peschiera Borromeo (MI).
<b>Sustainability:</b>	Environmental Management System complies with ISO 14001:2015.
<b>Materials:</b>	Tubing: PVC and silicone. All materials DEHP, BPA and Latex-free*. Additional details on request.
<b>Country of origin:</b>	Made in Italy.

\*DEHP and BPA are not intentionally added during the manufacturing process. Not formulated with natural rubber latex.

To obtain further information please contact your Compat® representative or go to [www.compat.com](http://www.compat.com)

®Unless stated otherwise, all trademarks are used with permission.

AXI00402 REV.1 16 MARCH 2021

# Compat®