

Compat ella[®]

Compat Ella[®] Spike Set

The Compat Ella[®] Spike Set is intended for administration of enteral nutrition, fluids or medication and is suitable for use with the Compat Ella[®] feeding pump.

Features at a Glance

The Compat Ella[®] Spike Set features an ENPlus nutrition connector, an ENFit[®] compatible 3-way medication port and an ENFit[®] distal end (patient side tube connection).


Reliable Feeding: The Compat Ella[®] Spike Set features patented fail-safe Auto-Anti-Free-Flow (AAFF) technology. AAFF is incorporated into the cassette and is always automatically activated to **ensure patient safety**.

Safer Connectors: The Compat Ella[®] Spike Set includes ENFit[®] connectors made according to ISO 80369-3 and an ENPlus connector made according to ISO 18250-3. ENFit[®] and ENPlus connectors are designed to be IV incompatible to **promote patient safety** by helping to avoid accidental misconnections.

Safe Materials: The Compat Ella[®] Spike Set is **sterile, DEHP, BPA and Latex-Free***.

Compatibility:

- Nutrition containers: ENPlus
- Enteral feeding tubes: ENFit[®]
- Syringes: ENFit[®]

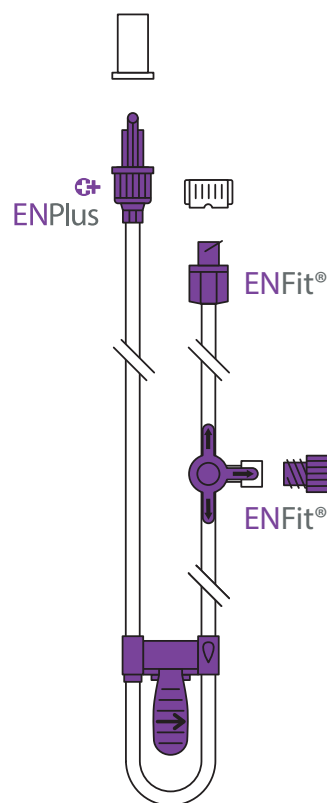


Ordering Information

Product Name	Product Code	Qty/Box
Compat Ella [®] Spike Set	12298203	30

Have you met Compat Ella[®]?

Compat Ella® Spike Set

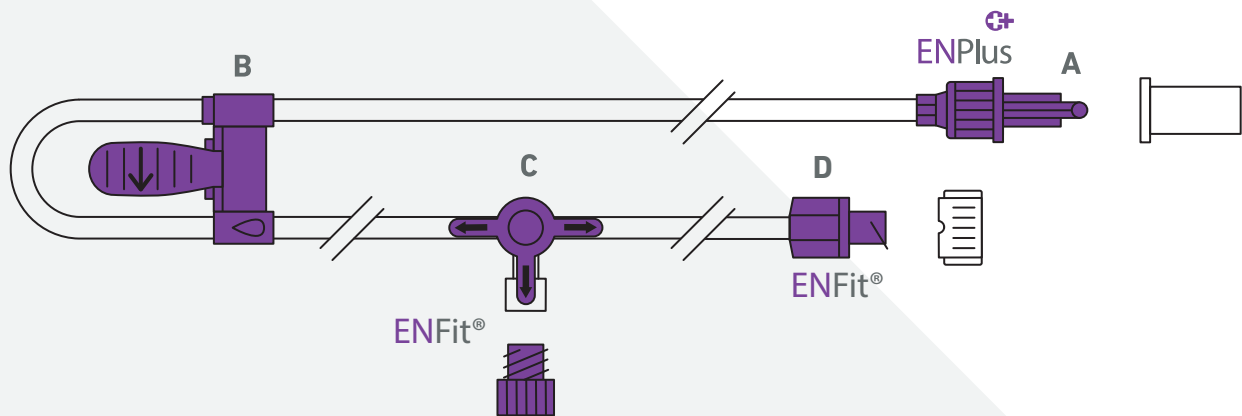

Components

- A ENPlus 'enteral only' connector (nutrition connection) with revolving lock and cap
- B Compat Ella® pump cassette with integrated and patented Auto-Anti-Free-Flow (AAFF) safety feature
- C ENFit® compatible 3-way medication port with cap
- D ENFit® distal end with cap (patient side tube connection)

Specifications

Sterile:	Sterilized with Ethylene Oxide (ETO).
Packing:	Individually packed in cases of 30 pieces (29.5 x 26 x 19 cm).
Storage:	Store at room temperature. Keep away from sunlight and humidity.
Shelf life:	3 years from production date. Expiry date on the package.
Duration of use:	24 hours. Single Use.
Length:	Distance from nutrition connection to patient side tube connection: 230 cm. Distance from nutrition connection to pump: 50 cm.
Diameters:	Nominal inner diameter: 3.0 mm. Nominal outer diameter: 4.15 mm.
ISO Standard:	ENFit® connector made according to ISO 80369-3 and ENPlus connector made according to ISO 18250-3.
CE Mark:	Complies with the European Medical Devices Directive 93/42/EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
Classification:	Class IIa.
Legal manufacturer:	Cedic Srl, Via Liberazione 63/9, IT-20068 Peschiera Borromeo (MI).
Sustainability:	Environmental Management System complies with ISO 14001:2015.
Materials:	Tubing: PVC and silicone. All materials DEHP, BPA and Latex-free*. Additional details on request.
Country of origin:	Made in Italy.

*DEHP and BPA are not intentionally added during the manufacturing process. Not formulated with natural rubber latex.

To obtain further information please contact your Compat® representative or go to www.compat.com

®Unless stated otherwise, all trademarks are used with permission.

AXI00405 REV.1 01 MARCH 2021

Compat®