

Compat[®] Flexibaggle[®]

Foldable Empty Enteral Reservoir 1.5 L

The Compat[®] Flexibaggle[®] 1.5 L foldable empty reservoir with ENPlus port is intended for the administration of enteral nutrition and/or fluids to patients. For use in conjunction with administration sets.

Features at a Glance


The Compat[®] Flexibaggle[®] 1.5 L foldable empty reservoir with ENPlus port has been designed to help improve patient care.

Ease of Use: The built-in hanger and large top-fill opening allow the bag to be easily filled with nutrition or water and its fold-ability saves on storage space enabling it to easily be at hand when needed.

Ambulatory Use: The bag resists a fall of up to 1.5 m and the easy-to-close cap design requires only a quarter turn to ensure a secure seal, allowing the product to also be used on-the-go with a low risk of leakage.

Patient Safety: Compat[®] Flexibaggle[®] 1.5 L includes an ENPlus connector made according to ISO 18250-3. The ENPlus connector is designed to be IV incompatible to promote patient safety by helping to avoid accidental misconnections. Compat[®] Flexibaggle[®] 1.5 L is made of Ethylene-vinyl Acetate (EVA), a soft and food compatible material. It is sterile, DEHP, BPA and Latex-Free*.

Eco-friendly: Compat[®] Flexibaggle[®] 1.5 L is foldable, reducing the environmental impact of shipping.

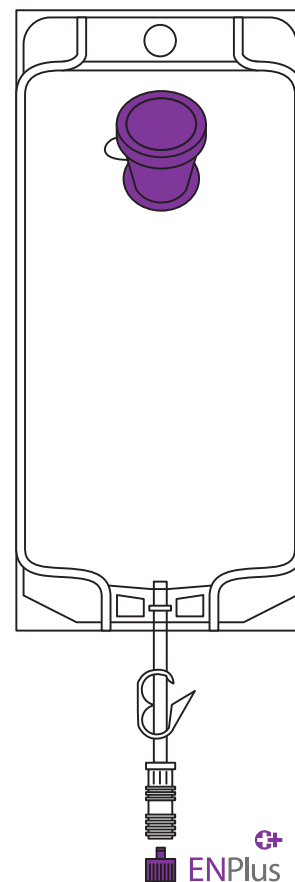


Ordering Information

Product Name	Product Code	Qty/Box
Compat [®] Flexibaggle [®] 1.5 L	12306162	30

Patient Care. Our Priority.

Foldable Empty Enteral Reservoir 1.5 L

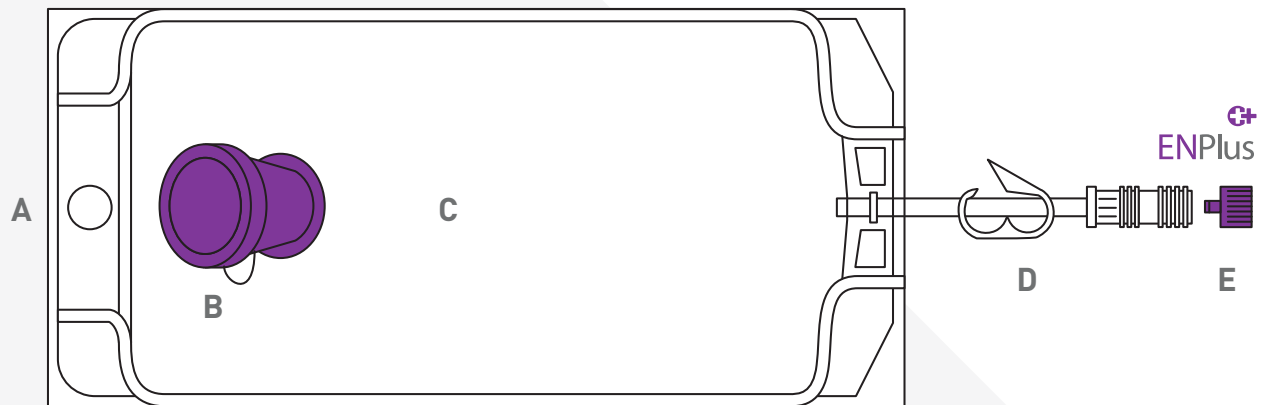

Components

- A Built-in hanger
- B Large top-fill opening with secure screw cap
- C 1.5 L container with graduations (at 50 ml and then every 100 ml) and tags for patient name and other important information (formula, volume, date and rate)
- D Pinch clamp
- E ENPlus 'enteral only' feeding port with cap

Specifications

- Wide-mouth port:** Nominal inner diameter: ± 36 mm.
- Bag size:** Max. nominal bag size: $\pm 326 \times 190$ mm.
- Sterile:** Sterilized with Ethylene Oxide (EO).
- Packing:** Individually packed in cases of 30 pieces (58.5 x 38.5 x 16 cm).
- Storage:** Store dry at room temperature and away from sunlight.
- Shelf life:** Up to 3 years from production date. Expiry date on the package.
- Duration of use:** 24 hours. Single Use.
- ISO Standard:** ENPlus connector made according to ISO 18250-3.
- CE Mark:** Complies with the European Medical Devices Directive 93/42/EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
- Classification:** Class Is.
- Legal manufacturer:** Cedic Srl, Via Liberazione 63/9, IT-20068 Peschiera Borromeo (MI).
- Sustainability:** Environmental management system complies with ISO 14001:2015.
- Materials:** Bag: EVA. Tubing: PVC. All materials DEHP, BPA and Latex-free*. Additional details on request.
- Country of origin:** Made in Italy.

*DEHP and BPA are not intentionally added during the manufacturing process. Not formulated with natural rubber latex.

To obtain further information please contact your Compat® representative or go to www.compat.com
*Unless stated otherwise, all trademarks are used with permission.

AXI00032 REV.1 08 MARCH 2021

Compat®