

Compat[®] Gravity

Compat[®] Gravity Universal Set

The Compat[®] Gravity Universal Set is intended for administration of enteral nutrition, fluids or medication and is suitable for gravity use.

Features at a Glance


The Compat[®] Gravity Universal Set features an ENPlus nutrition connector, a universal screw cap with air filter, a drip chamber, a roller clamp, an ENFit[®] Y medication port and an ENFit[®] distal end. A separate bottle hanger is also included.

Patient Safety: Compat[®] Gravity Universal Set includes ENFit[®] connectors made according to ISO 80369-3 and ENPlus connectors made according to ISO 18250-3. The ENFit[®] and ENPlus connectors are designed to be IV incompatible to **promote patient safety** by helping to avoid accidental misconnections.

Safe Materials: Compat[®] Gravity Universal Set is **sterile, DEHP, BPA and Latex-Free***.

Compatibility:

- Nutrition containers: ENPlus, glass and plastic wide-neck containers, crown cork glass bottles
- Enteral feeding tubes: ENFit[®]
- Syringes: ENFit[®]

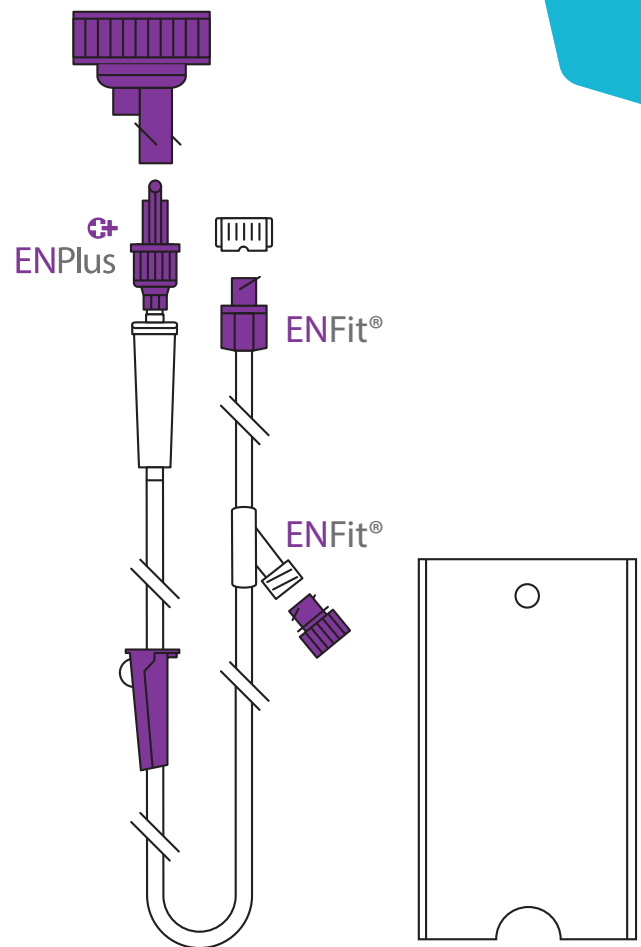


Ordering Information

Product Name	Product Code	Qty/Box
Compat [®] Gravity Universal Set	12298230	30

Patient Care. Our Priority.

Compat[®] Gravity Universal Set

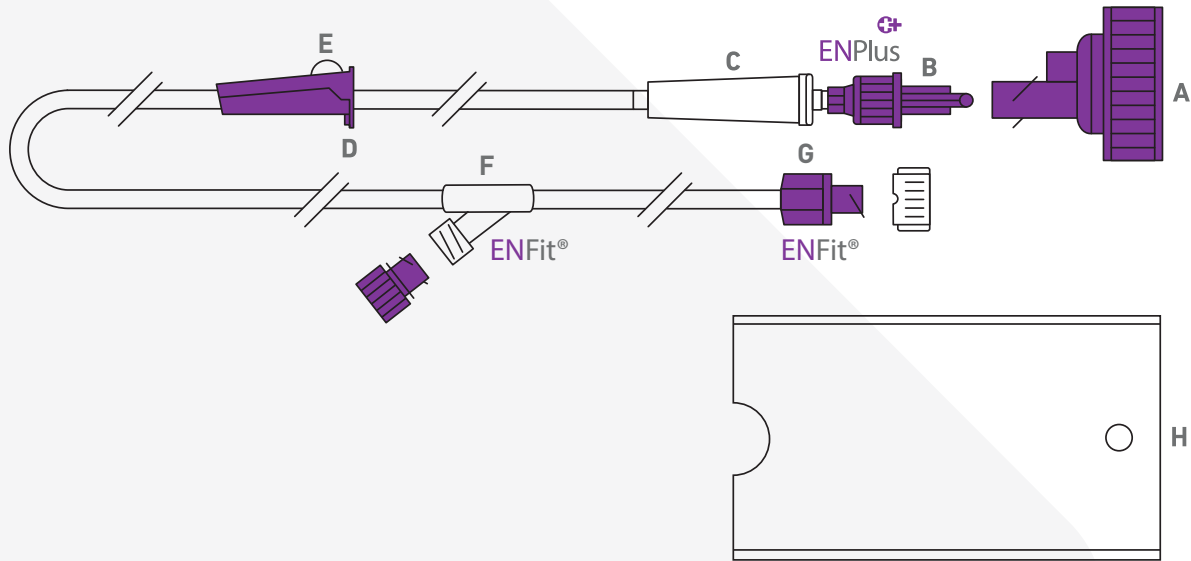

Components

- A** Universal screw cap with air filter and ENPlus port compatible with glass and plastic Ø40 mm wide-neck containers (DIN 168-1) or crown cork glass bottles (EN 14634 and EN 14635)
- B** ENPlus 'enteral only' proximal end (nutrition connection) with revolving lock
- C** Drip chamber
- D** Built-in hanger to help keep the administration set distal end away from dirt
- E** Roller clamp to shut off or adjust the flow rate
- F** ENFit[®] Y medication port with cap
- G** ENFit[®] distal end with cap (patient side tube connection)
- H** 500 ml bottle holder LDPE (Low Density Polyethylene)

Specifications

- Sterile:** Sterilized with Ethylene Oxide (EO).
- Packing:** Individually packed in cases of 30 pieces (39 x 39 x 13 cm).
- Storage:** Store at room temperature. Keep away from sunlight and humidity.
- Shelf life:** 3 years from production date. Expiry date on the package.
- Duration of use:** 24 hours. Single Use.
- Tube Length:** Distance from nutrition connection to patient side tube connection: 180 cm.
- Tube diameters:** Nominal inner diameter: 4.0 mm. Nominal outer diameter: 5.4 mm.
- ISO Standard :** ENFit[®] connector made according to ISO 80369-3 and ENPlus connector made according to ISO 18250-3.
- CE Mark:** Complies with the European Medical Devices Directive 93/42/EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
- Classification:** Class Is.
- Legal manufacturer:** Cedic Srl, Via Liberazione 63/9, IT-20068 Peschiera Borromeo (MI).
- Sustainability:** Environmental management system complies with ISO 14001:2015.
- Materials:** Tubing: PVC. All materials DEHP, BPA and Latex-free*. Additional details on request.
- Country of origin:** Made in Italy.

*DEHP and BPA are not intentionally added during the manufacturing process. Not formulated with natural rubber latex.

To obtain further information please contact your Compat[®] representative or go to www.compat.com
®Unless stated otherwise, all trademarks are used with permission.

AXI00036 REV.1 08 MARCH 2021

Compat[®]