

# Compat<sup>®</sup> Modum

## Closed System Gastric Residuals Aspiration and Measurement Accessory

Compat<sup>®</sup> Modum is intended for use as an accessory to enable closed system gastric residuals aspiration and their possible re-introduction into the stomach of the patient, when used in conjunction with appropriate tubes, syringes and collection bags. Also intended for use for enteral feeding and medication delivery.

### Features at a Glance

Compat<sup>®</sup> Modum is an accessory placed between the feeding tube and the administration set to **ease gastric fluid aspiration and GRV\* measurement**.


#### Ease of Use:

- Compat<sup>®</sup> Modum enables **closed system** gastric residuals aspiration into a drainage collection bag **reducing exposure to gastric fluids and their smell and lowering the risk of contamination**.
- It features two 3-way connectors which allow to **select the flow path** to easily switch from gastric emptying to enteral feeding or medication administration.
- It is **compatible** with most enteral feeding tubes, medication dispensers and drainage/suction devices\*\*.

**Safe Materials:** Compat<sup>®</sup> Modum is **sterile, DEHP, BPA and Latex-Free\*\*\***.

#### Compatibility:

- Administration sets: ENFit<sup>®</sup>
- Enteral feeding tubes: ENFit<sup>®</sup> and ENLock
- Syringes: ENFit<sup>®</sup>, ENLock and Luer
- Drainage/suction devices: Funnel\*\*

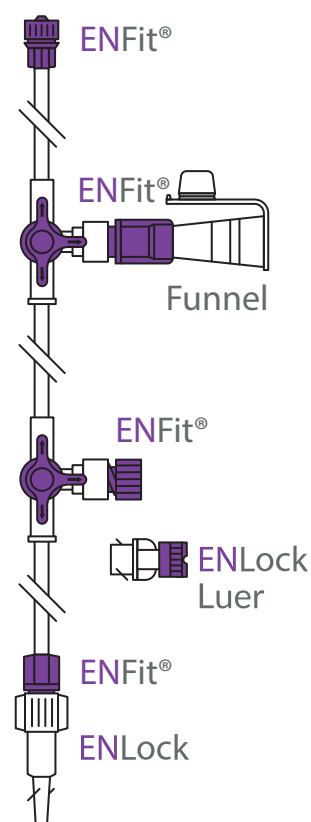


### Ordering Information

Product Name	Product Code	Qty/Box
Compat <sup>®</sup> Modum	12340186	50

*Patient Care. Our Priority.*

# Closed System Gastric Residuals Aspiration and Measurement Accessory

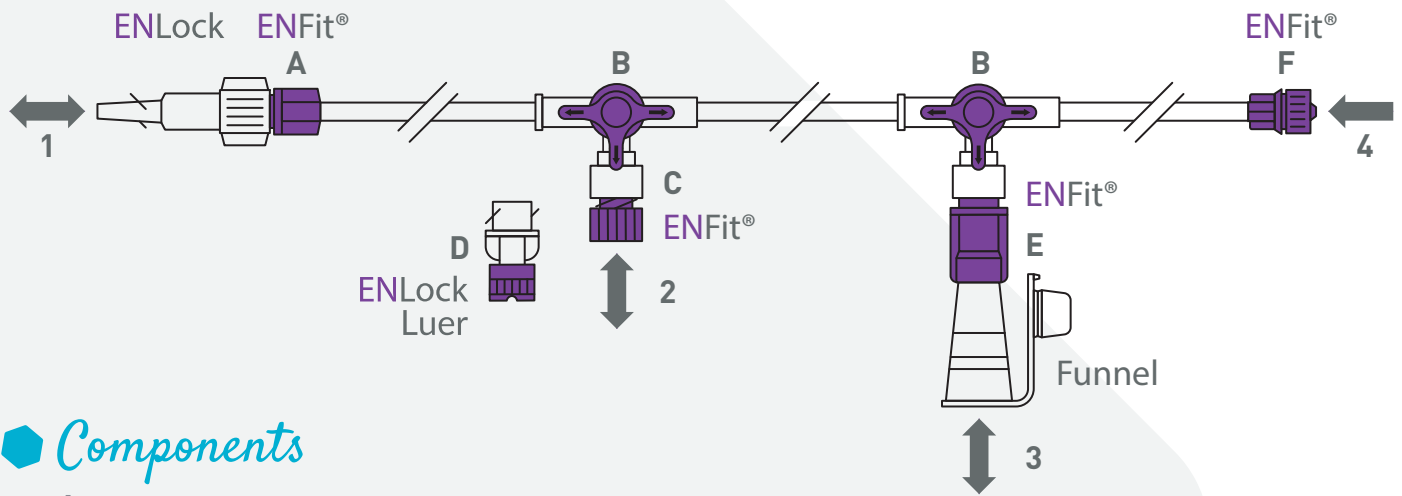

## Components

- A Female ENFit® distal end with removable ENLock adapter
- B ENFit® compatible 3-way connectors to select the flow path
- C Male ENFit® compatible drainage syringe port/medication port with cap
- D ENFit® adapter with cap compatible with Luer, Oral and ENLock medication dispensers/drainage syringes
- E Funnel connector with cap designed to be compatible with most drainage/suction devices\*\*
- F Male ENFit® proximal connector for connection to administration set

## Flow Paths

- 1 To or from feeding tube (patient side connection)
- 2 To or from drainage syringe/From medication dispenser
- 3 To or from drainage/suction devices
- 4 From administration set (nutrition side connection)

## Specifications

<b>Sterile:</b>	Sterilized with Ethylene Oxide (EO).
<b>Length:</b>	12 cm.
<b>Packing:</b>	Individually packed in cases of 50 pieces (295x260x190mm).
<b>Storage:</b>	Store dry at room temperature and away from sunlight.
<b>Shelf life:</b>	Up to 3 years from production date. Expiry date on the package.
<b>Duration of use:</b>	Up to 7 days. Single Use.
<b>ISO Standard:</b>	ENFit® connector made according to ISO 80369-3.
<b>CE Mark:</b>	Complies with the European Medical Devices Directive 93/42/EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
<b>Classification:</b>	Class IIa.
<b>Legal manufacturer:</b>	Cedic Srl, Via Liberazione 63/9, IT-20068 Peschiera Borromeo (MI).
<b>Sustainability:</b>	Environmental management system complies with ISO 14001:2015.
<b>Materials:</b>	Tubing: PVC. All materials DEHP, BPA and Latex-free***.
<b>Country of origin:</b>	Made in Italy.

\*GRV: Gastric Residual Volume.

\*\*Some incompatibilities may exist.

\*\*\*DEHP and BPA are not intentionally added during the manufacturing process. Not formulated with natural rubber latex.

To obtain further information please contact your Compat® representative or go to [www.compat.com](http://www.compat.com)

®Unless stated otherwise, all trademarks are used with permission.

AXI00037 REV.3 16 MARCH 2021

# Compat®