

Compat[®] Buddy

Full Range of Low-profile Gastrostomy Feeding Tubes

The Compat[®] Buddy low-profile gastrostomy feeding tube (button) is intended for use in patients who require intermittent or continuous enteral administration of nutrition, fluids or medication at the gastric site. Also suitable for gastric decompression.

Features at a Glance


The Compat[®] Buddy button is available in a variety of sizes (12/14/16/18/20/24FR) and lengths (0.8 to 6.5 cm) to address multiple patient and clinician needs.

Patient Comfort:

- Soft biocompatible silicone tube and high-resistance silicone balloon to promote maximum service life and limit tube replacement needs.
- With recessed distal end to be gentle on the stomach wall.
- Discreet and comfortable solution to offer greater freedom of movement to the patient.

Ease of Use: The Compat[®] Buddy button is available as a convenient kit for easy placement. The kit is delivered in individual, sterile packaging with two detachable labels to improve traceability.

Safe Materials: DEHP, BPA and Latex-free*.

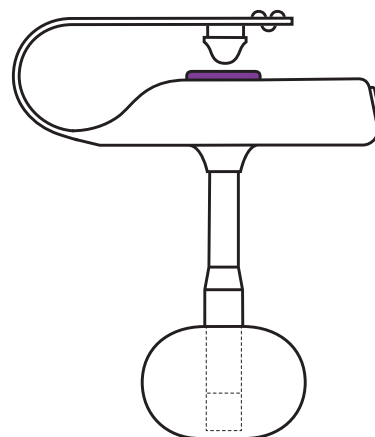


Ordering Information

Stoma tract length (cm)	French size					
	12FR	14FR	16FR	18FR	20FR	24FR
0.8	12466557	12466601	12466635	12466759	12466828	12466844
1.0	12466535	12466577	12466642	12466795	12466835	12466864
1.2	12466544	12466600	12466636	12466784	12466841	12466847
1.5	12460321	12460325	12466761	12466785	12466836	12466870
1.7	12466556	12460330	12460320	12466798	12466815	12466865
2.0	12466534	12460327	12460341	12466823	12466842	12466846
2.3	12466543	12466568	12466745	12466824	12466817	12466856
2.5	12466555	12460342	12460352	12466825	12466837	12466871
2.7	12466554	12466576	12466719	12466805	12466818	12466857
3.0	12466561	12466589	12460347	12466830	12466838	12466869
3.5	12466560	12460331	12466752	12466806	12466839	12466868
4.0	12466551	12466620	12466753	12466826	12466863	12466859
4.5	12466528	12466598	12466770	12466807	12466862	12466890
5.0	12466527	12466621	12466762	12466808	12466852	12466892
5.5		12466616		12466811		
6.0		12466632		12466809		
6.5		12466626				

Patient Care. Our priority.

Full Range of Low-profile Gastrostomy Feeding Tubes

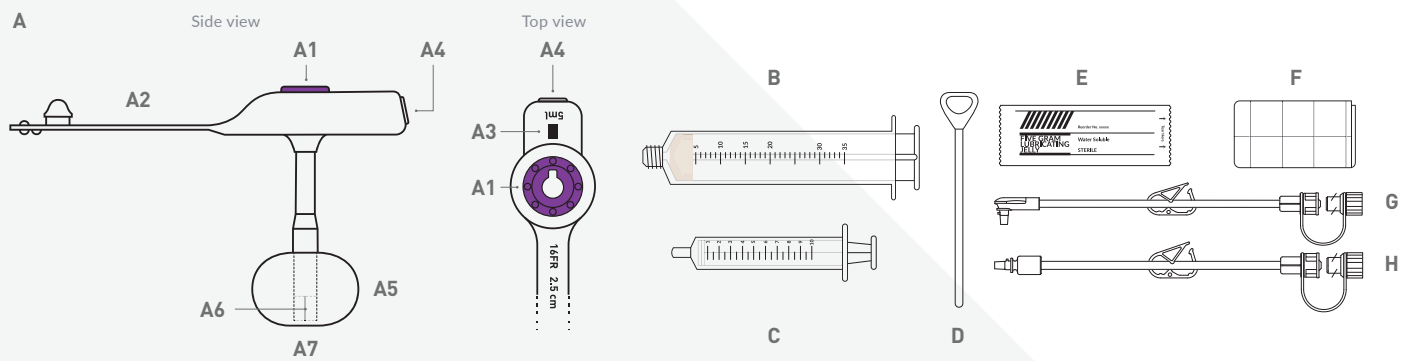

Components

- A** Compat® Buddy low-profile gastrostomy feeding tube
A1 Purple-colored feeding port with anti-reflux valve. Compatible with straight and right-angle extension sets
A2 Cap with indication of the French size and length
A3 Position mark for the extension set
A4 Inflation valve with Luer connector and indication of balloon volume, in ml

Button size and balloon volume						
Size (FR)	12	14	16	18	20	24
Recommended inflation volume (ml)	3.0	4.0	5.0	5.0	8.0	8.0
Max. inflation volume (ml)	5.0	6.0	7.0	7.0	12.0	12.0

- A5** High-resistance silicone balloon
A6 Radio-opaque stripe for placement and position confirmation
A7 Recessed distal tip
B ENFit® 35 ml syringe compatible with ENFit® extension sets
C Luer syringe for balloon inflation
D Stylet
E Water-soluble lubricant
F Gauze pads (2 pieces)
G 30 cm single ENFit® extension set with right-angle connector
H 30 cm single ENFit® extension set with straight connector

Specifications

- Sterile:** Sterilized with Ethylene Oxide (EO).
Packaging: Individually packed.
Storage: Store in the original packaging at room temperature and in a dry environment. Avoid direct sunlight and heat radiation.
Shelf life: 3 years from production date. Expiry date on the package.
Duration of use: Indicated for long-term use. For single use.
Size and length: From 12 to 24FR and from 0.8 to 6.5 cm. See table on page 1.
CE Mark: Complies with the European Medical Devices Directive 93/42/EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
Classification: Class IIb.
Legal manufacturer: Degania Silicone Ltd.
Materials: Tube: Silicone. All materials are DEHP, BPA and Latex-free*. Additional details on request.

*DEHP and BPA are not intentionally added during the manufacturing process. Not formulated with natural rubber latex.

To obtain further information please contact your Compat® representative or go to www.compat.com
 *Unless otherwise stated, all trademarks are used with permission.