

# Compat<sup>®</sup> DualPort

## Single Lumen Dual Port Nasogastric Feeding and Drainage Tube

Compat<sup>®</sup> DualPort is a single lumen, dual port nasogastric tube for: 1) the administration of enteral nutrition, fluids and medication via the ENFit<sup>®</sup> enteral-specific feeding port and 2) gastric drainage/intermittent aspiration of stomach contents via the funnel type drainage port.

### Features at a Glance

**Designed for the ICU:** Compat<sup>®</sup> DualPort features a funnel type as well as an ENFit<sup>®</sup> connector to help simplify gastric drainage and administration of nutrition, fluids and medication by the use of one tube for both operations.

- Allows both gastric feeding and drainage/aspiration via a single tube.
- Prevents the need for a second tube placement when the patient transitions from requiring drainage/aspiration to enteral feeding.
- Designed to resist collapsing when checking gastric residuals.


Clamps help prevent fluid leakage and ease gastric fluid aspiration.

**Patient Comfort:** Compat<sup>®</sup> DualPort is made of soft radio-opaque polyurethane with a rounded distal end to help promote patient comfort.

**Aid in Placement:** A pre-lubricated and open “flow-through” stylet is included to aid in placement:

- Allows injecting water and radio-opaque fluids without stylet removal.
- May help lower the risk of tube displacement on stylet removal.

Compat<sup>®</sup> DualPort is compatible with all ENFit<sup>®</sup> administration sets and syringes as well as most drainage/suction devices\*.

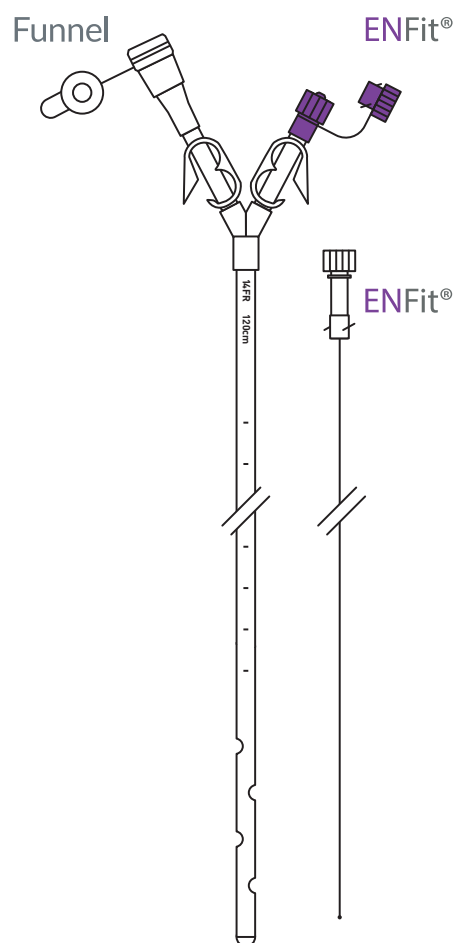


### Ordering Information

Product Name	Product Code	Qty/Box
Compat <sup>®</sup> DualPort 14FR 120 CM S	12391311	20

*Patient Care. Our Priority.*

# Single Lumen Dual Port Nasogastric Feeding and Drainage Tube

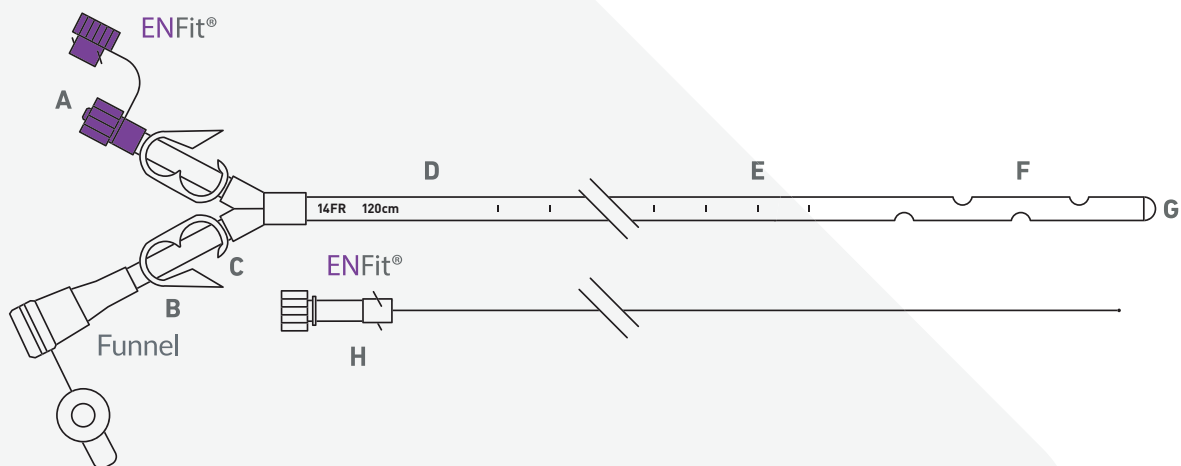

## Components

- A ENFit® twist-lockable feeding connector with cap
- B Funnel connector with cap designed to be compatible with most drainage/suction devices\*
- C Closing clamps
- D Product identification, printed French (FR) size and length on the shaft
- E Radio-opaque polyurethane tube shaft with placement depth markers every 5 cm
- F 4 lateral alternate openings
- G Rounded distal end
- H Pre-lubricated and open “flow-through” stylet

## Specifications

- Sterile:** Sterilized with Ethylene Oxide (EO).
- Packing:** Individually packed in cases of 20 pieces (29.5 × 25 × 17 cm).
- Storage:** Store in the original packaging at room temperature and in a dry environment. Avoid direct sunlight and heat radiation.
- Shelf life:** 3 years from production date. Expiry date on the package.
- Duration of use:** Up to 29 days. Single Use.
- Length:** 120 cm.
- Diameter:** 14FR.
- ISO Standard:** ENFit® connector made according to ISO 80369-3.
- CE Mark:** Complies with the European Medical Devices Directive 93/42/EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
- Classification:** Class IIa.
- Legal manufacturer:** Cedic Srl, Via Liberazione 63/9, IT-20068 Peschiera Borromeo (MI).
- Sustainability:** Environmental management system complies with ISO 14001:2015.
- Materials:** Tubing: Polyurethane. ENFit® connector: Copolyester. Funnel connector: Polyurethane. DEHP, BPA and Latex-free\*\*. Additional details on request.
- Country of origin:** Made in Italy.

\*Some incompatibilities may exist.  
\*\*DEHP and BPA are not intentionally added during the manufacturing process. Not formulated with natural rubber latex.  
To obtain further information please contact your Compat® representative or go to [www.compat.com](http://www.compat.com)  
\*Unless stated otherwise, all trademarks are used with permission.