

Compat[®] Buddy

Full Range of Low-profile Gastrostomy Feeding Tube Extension Sets

The Compat[®] Buddy Low-profile Gastrostomy Feeding Tube Extension Sets are intended for the administration of enteral nutrition, fluids and/or medication. Indicated for patients who require nutritional support and have a low-profile gastrostomy feeding tube (button). Also suitable for passive gastric decompression.

Features at a Glance

Multiple Sizes and Features:


- Available with right-angle and straight connector.
- 5, 30 and 60 cm to address multiple patient and clinician needs.
- Dual port to facilitate medication administration (optional).
- Screw cap and closing clamp (optional) to help prevent fluid leakage and ease handling.

Compatibility:

- Compat[®] Buddy Extension Sets are compatible with Compat[®] Buddy button and with most other buttons including Mic-Key[®], Aquarius[™], MiniONE[®], Freka[®] Belly, Tummy Button[®] and danuButton[®].
- Only compatible with administration sets and syringes with ENFit[®] connector to promote patient safety.

Patient Safety: ENFit[®] connectors are made according to ISO 80369-3. They are designed to be IV incompatible to promote patient safety by helping to avoid accidental misconnections.

Safe Materials: DEHP, BPA and Latex-Free*.

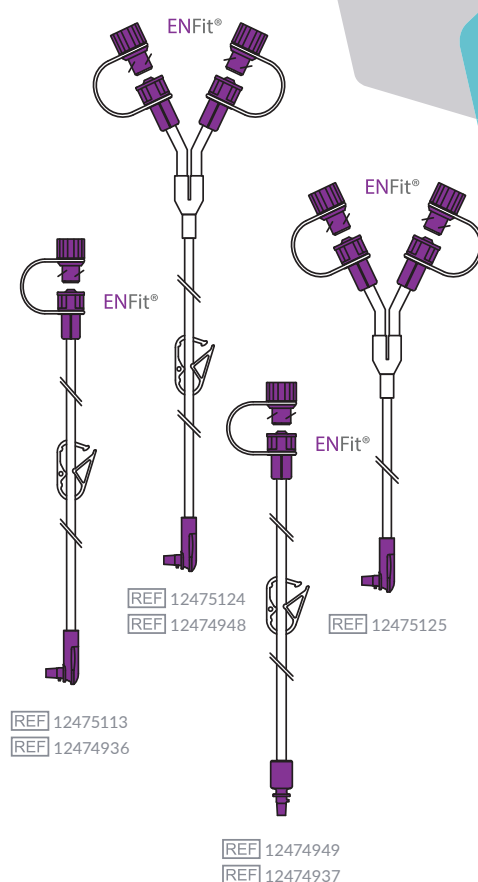


Ordering Information

	Connector	ENFit [®] Port	Length (cm)	Product Code	Qty/box
Compat [®] Buddy Low-profile Gastrostomy Feeding Tube Extension Set	Right-angle	Dual	5	12475125	10
	Right-angle	Dual	30	12475124	
	Right-angle	Dual	60	12474948	
	Right-angle	Single	30	12475113	
	Right-angle	Single	60	12474936	
	Straight	Single	30	12474949	
	Straight	Single	60	12474937	

Patient Care. Our Priority.

Full Range of Low-profile Gastrostomy Feeding Tube Extension Sets

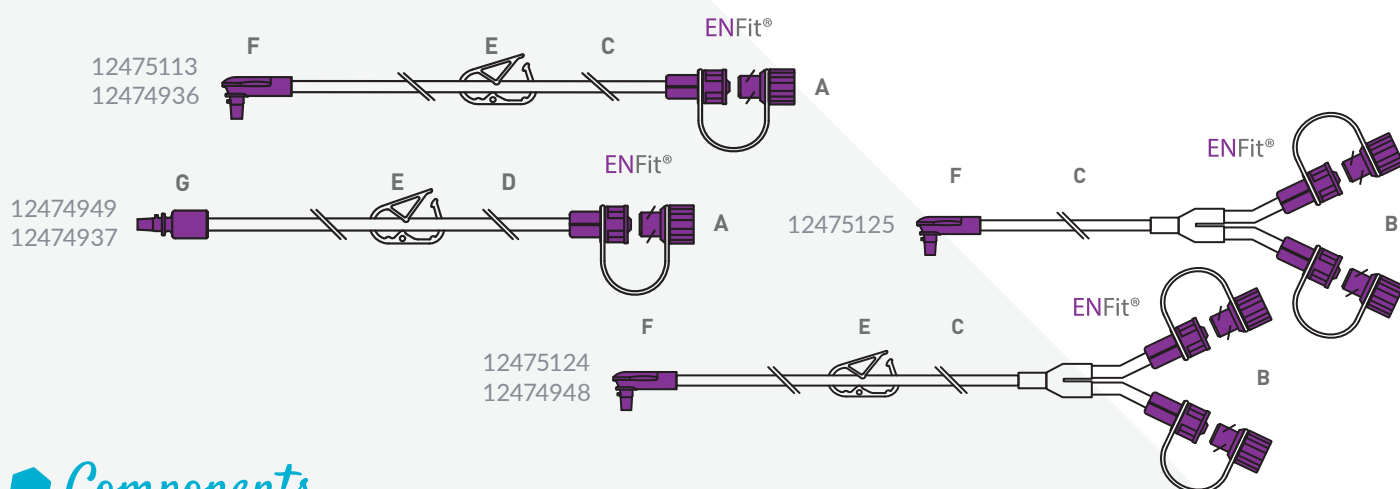

Components

- A** ENFit® twist-lockable feeding port with cap and tether
- B** Dual ENFit® twist-lockable feeding port with cap and tether
- C** Standard PVC tubing outer diameter 4.15 mm
- D** Large bore PVC tubing outer diameter 5.4 mm
- E** Closing clamp
- F** Right-angle connector with alignment marker (patient side tube connection)
- G** Straight connector with alignment marker (patient side tube connection)

Specifications

- Sterile:** Sterilized with Ethylene Oxide (EO).
- Packing:** Individually packed in a box of:
10 pieces per case (19 × 19 × 10 cm).
24 cases in a shipping box (60 × 40 × 42 cm).
- Storage:** Store dry at room temperature and away from sunlight.
- Shelf life:** 3 years from production date. Expiry date on the package.
- Duration of use:** Up to 7 days and 5 usages per day.
- Tube length:** 5, 30 or 60 cm.
- Outer tube diameter:** 4.15 mm (REF with right-angle connector) or 5.4 mm (REF with straight connector).
- ISO Standard:** ENFit® connector made according to ISO 80369-3.
- CE mark:** Complies with the European Medical Devices Directive 93/42/EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
- Classification:** Class IIa.
- Legal manufacturer:** Cedec Srl, Via Liberazione 63/9, IT-20068 Peschiera Borromeo (MI).
- Sustainability:** Environmental management system complies with ISO 14001:2015.
- Materials:** Tubing: PVC. ENFit® and button connectors: Copolyester. All materials DEHP, BPA and Latex-free*. Additional details on request.
- Country of origin:** Made in Italy.

*DEHP and BPA are not intentionally added during the manufacturing process. Not formulated with natural rubber latex.
ENFit® and Compat® trademarks are used with permission. Mic-Key® is a registered trademark of Halyard Health, Inc. or its affiliates.
Freka® is a registered trademark of Fresenius Kabi AG. MiniONE® is a registered trademark of Applied Medical Technology,
Aquarius™ is a registered trademark of Degania Medical, Inc. Tummy Button® is a registered trademark of MediLime B.V.
danuButton® is a registered trademark of danumed Medizintechnik GmbH.

To obtain further information please contact your Compat® representative or go to www.compat.com

AXI00469 REV.5 17 JUNE 2024

Compat®