

Compat[®] Modum

Closed System Gastric Residuals Aspiration Accessory

Compat[®] Modum is intended for use as an accessory to enable closed system gastric residuals aspiration and their possible re-introduction into the stomach of the patient, when used in conjunction with appropriate tubes, syringes and collection bags. Also intended for use for enteral feeding and medication delivery.

Features at a Glance

Compat[®] Modum is an accessory placed between the feeding tube and the administration set to **ease gastric fluid aspiration for GRV* measurement**.


Ease of Use:

- Compat[®] Modum enables **closed system** gastric residuals aspiration into a drainage collection bag **reducing exposure to gastric fluids and their smell and lowering the risk of contamination**.
- It features two 3-way connectors which allow to **select the flow path** to easily switch from gastric emptying to enteral feeding or medication administration.
- It is **compatible** with most enteral feeding tubes, medication dispensers and drainage/suction devices**.

Safe Materials: Compat[®] Modum is **sterile, DEHP, BPA and Latex-Free*****.

Compatibility:

- Administration sets: ENFit[®]
- Enteral feeding tubes: ENFit[®] and ENLock
- Syringes: ENFit[®], ENLock and Luer
- Drainage/suction devices: Funnel**

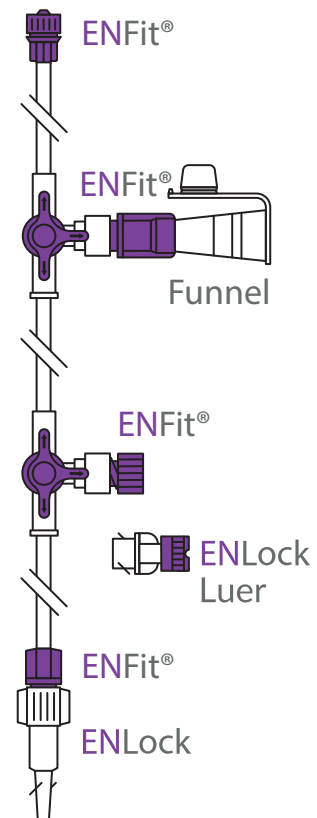


Ordering Information

Product Name	Product Code	Qty/Box
Compat [®] Modum	12340186	50

Patient Care. Our Priority.

Closed System Gastric Residuals Aspiration Accessory

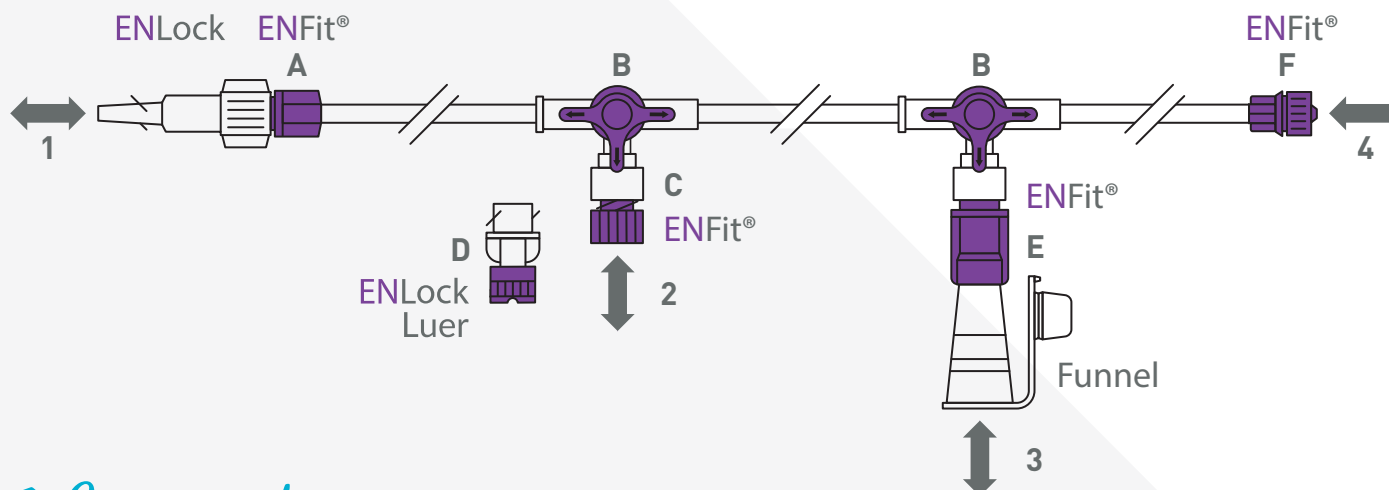

Components

- A Female ENFit® distal end with removable ENLock adapter
- B ENFit® 3-way connectors to select the flow path
- C Male ENFit® drainage syringe port/medication port with cap
- D ENFit® adapter with cap compatible with Luer, Oral and ENLock medication dispensers/drainage syringes
- E Funnel connector with cap designed to be compatible with most drainage/suction devices**
- F Male ENFit® proximal connector for connection to administration set

Flow Paths

- 1 To or from feeding tube (patient side connection)
- 2 To or from drainage syringe/from medication dispenser
- 3 To or from drainage/suction devices
- 4 From administration set (nutrition side connection)

Specifications

Sterile:	Sterilized with Ethylene Oxide (EO).
Length:	12 cm.
Packing:	Individually packed in cases of 50 pieces (295 × 260 × 190 mm).
Storage:	Store dry at room temperature and away from sunlight.
Shelf life:	Up to 3 years from production date. Expiry date on the package.
Duration of use:	Up to 7 days. Single Use.
ISO Standard:	ENFit® connector made according to ISO 80369-3.
CE Mark:	Complies with the European Medical Devices Regulation 2017/745. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
Classification:	Class IIa.
Legal manufacturer:	Cedic Srl, Via Liberazione 63/9, IT-20068 Peschiera Borromeo (MI).
Sustainability:	Environmental management system complies with ISO 14001.
Materials:	Tubing: PVC. All materials DEHP, BPA and Latex-free***.
Country of origin:	Made in Italy.

*GRV: Gastric Residual Volume.

**Some incompatibilities may exist.

***DEHP and BPA are not intentionally added during the manufacturing process. Not formulated with natural rubber latex.

This product datasheet does not replace the original Instructions for Use which must be read before use. To obtain further information please contact your Compat® representative or go to www.compat.com

®Unless stated otherwise, all trademarks are used with permission.

AXI00037 REV.4 03 OCTOBER 2025

Compat®