

Compat ella[®]

Compat Ella[®] Combiset 1.5 L

The Compat Ella[®] Combiset 1.5 L is intended for administration of enteral nutrition, fluids or medication and is suitable for use with the Compat Ella[®] feeding pump.

Features at a Glance

The Compat Ella[®] Combiset 1.5 L features an integrated 1.5 L Flexibaggle[®] reservoir, an ENFit[®] 3-way medication port and an ENFit[®] distal end (patient side tube connection).

Reliable Feeding: The Compat Ella[®] Combiset 1.5 L features patented fail-safe Auto-Anti-Free-Flow (AAFF) technology. AAFF is incorporated into the cassette and is always automatically activated to ensure patient safety.


Safer Connectors: The Compat Ella[®] Combiset 1.5 L includes ENFit[®] connectors made according to ISO 80369-3. The ENFit[®] connector is designed to be IV incompatible to promote patient safety by helping to avoid accidental misconnections.

Convenience: The integrated 1.5 L reservoir offers a convenient open system with a built-in hanger and screw cap. The large top-fill opening allows the bag to be refilled with nutrition or water without stopping the feeding program, or removing the reservoir from the hanger.

Safe Materials: The 1.5 L reservoir is made of Ethylene-vinyl Acetate (EVA), a soft and food compatible material. The Compat Ella[®] Combiset 1.5 L is sterile, DEHP, BPA and Latex-Free*.

Compatibility:

- Enteral feeding tubes: ENFit[®]
- Syringes: ENFit[®]

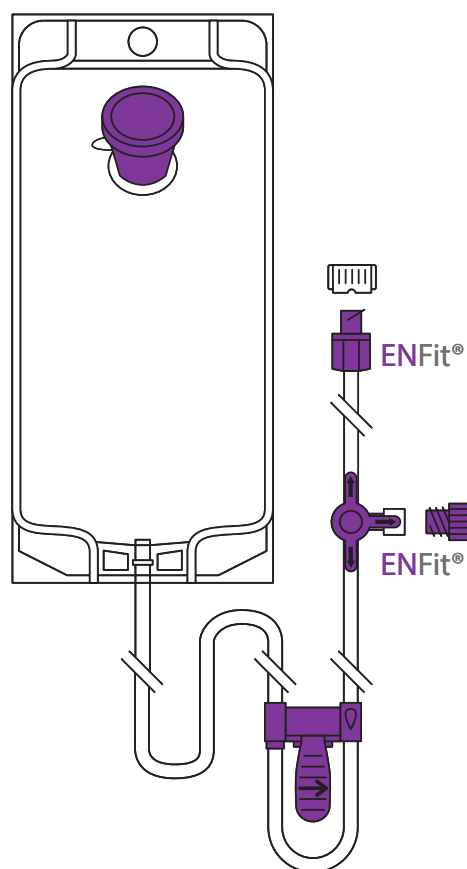


Ordering Information

Product Name	Product Code	Qty/Box
Compat Ella [®] Combiset 1.5 L	12298202	30

Have you met Compat Ella[®]?

Compat Ella® Combiset 1.5 L

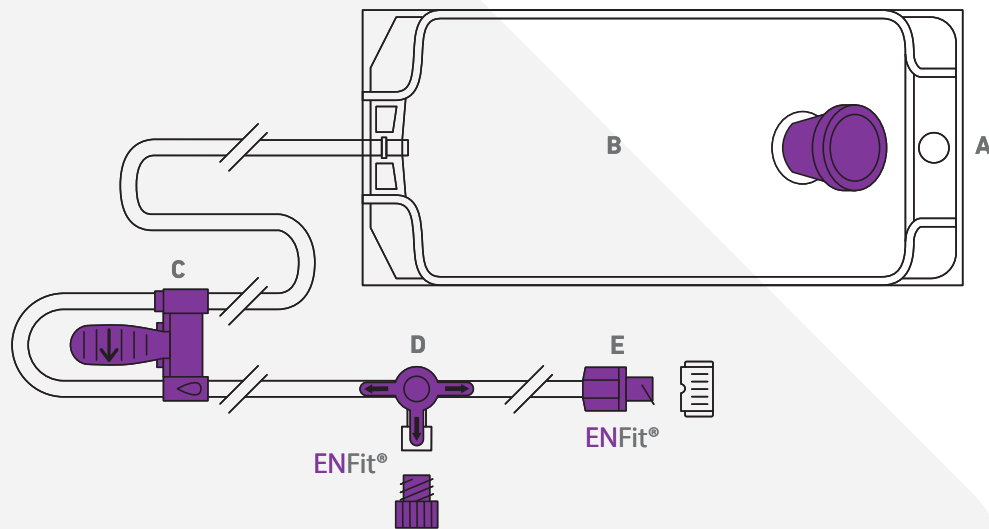

Components

- A** Large top-fill opening with built-in hanger and screw cap
- B** Flexibaggle® 1.5 L container made of EVA with graduations (at 50 ml and every 100 ml) and tags for patient name and other important information (formula, volume, date and rate)
- C** Compat Ella® pump cassette with integrated and patented Auto-Anti-Free-Flow (AAFF) safety feature
- D** ENFit® compatible 3-way medication port with cap
- E** ENFit® distal end with cap (patient side tube connection)

Specifications

Sterile:	Sterilized with Ethylene Oxide (EO).
Packing:	Individually packed in cases of 30 pieces (58.5 x 38.5 x 16 cm).
Storage:	Store at room temperature. Keep away from sunlight and humidity.
Shelf life:	3 years from production date. Expiry date on the package.
Duration of use:	24 hours. Single Use.
Length:	Distance from bag to patient side tube connection: 230 cm. Distance from Flexibaggle® to pump: 50 cm.
Diameters:	Nominal inner diameter: 3.0 mm. Nominal outer diameter: 4.15 mm.
ISO Standard:	ENFit® connector made according to ISO 80369-3.
CE Mark:	Complies with the European Medical Devices Directive 93/42/EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
Classification:	Class IIa.
Legal manufacturer:	Cedic Srl, Via Liberazione 63/9, IT-20068 Peschiera Borromeo (MI).
Sustainability:	Environmental Management System complies with ISO 14001.
Materials:	Tubing: PVC and silicone. All materials DEHP, BPA and Latex-free*. Additional details on request.
Country of origin:	Made in Egypt.

*DEHP and BPA are not intentionally added during the manufacturing process. Not formulated with natural rubber latex.

To obtain further information please contact your Compat® representative or go to www.compat.com

®Unless stated otherwise, all trademarks are used with permission.

AXI00402 REV.2 02 FEBRUARY 2026

Compat®