

Compat[®] Flexibaggle[®]

Foldable Empty Enteral Reservoir 1.5 L

The Compat[®] Flexibaggle[®] 1.5 L foldable empty reservoir with ENPlus port is intended for the administration of enteral nutrition and/or fluids to patients. For use in conjunction with administration sets.

Features at a Glance


The Compat[®] Flexibaggle[®] 1.5 L foldable empty reservoir with ENPlus port has been designed to help improve patient care.

Ease of Use: The built-in hanger and large top-fill opening allow the bag to be easily filled with nutrition or water and its foldability saves on storage space, enabling it to easily be at hand when needed.

Ambulatory Use: The bag resists a fall of up to 1.5 m and the easy-to-close cap design requires only a quarter turn to ensure a secure seal, allowing the product to also be used on-the-go with a low risk of leakage.

Patient Safety: Compat[®] Flexibaggle[®] 1.5 L includes an ENPlus connector made according to ISO 18250-3. The ENPlus connector is designed to be IV incompatible to promote patient safety by helping to avoid accidental misconnections. Compat[®] Flexibaggle[®] 1.5 L is made of Ethylene-vinyl Acetate (EVA), a soft and food compatible material. It is sterile, DEHP, BPA and Latex-Free*.

Eco-friendly: Compat[®] Flexibaggle[®] 1.5 L is foldable, reducing the environmental impact of shipping.

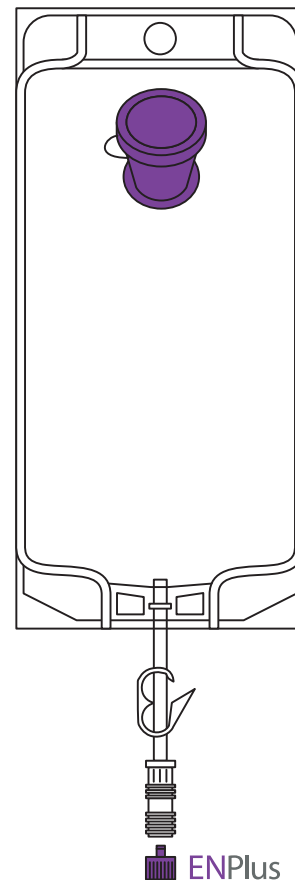


Ordering Information

Product Name	Product Code	Qty/Box
Compat [®] Flexibaggle [®] 1.5 L	12306162	30

Patient Care. Our Priority.

Foldable Empty Enteral Reservoir 1.5 L

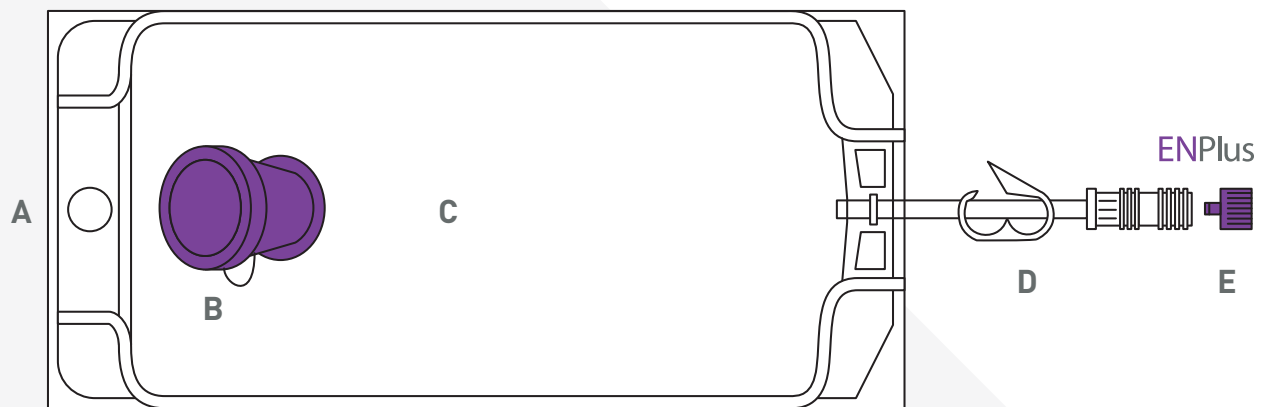

Components

- A Built-in hanger
- B Large top-fill opening with secure screw cap
- C 1.5 L container with graduations (at 50 ml and then every 100 ml) and tags for patient name and other important information (formula, volume, date and rate)
- D Pinch clamp
- E ENPlus 'enteral only' feeding port with cap

Specifications

- Wide-mouth port:** Nominal inner diameter: ± 36 mm.
- Bag size:** Max. nominal bag size: $\pm 326 \times 190$ mm.
- Sterile:** Sterilized with Ethylene Oxide (EO).
- Packing:** Individually packed in cases of 30 pieces (58.5 x 38.5 x 16 cm).
- Storage:** Store dry at room temperature and away from sunlight.
- Shelf life:** Up to 3 years from production date. Expiry date on the package.
- Duration of use:** 24 hours. Single Use.
- ISO Standard:** ENPlus connector made according to ISO 18250-3.
- CE Mark:** Complies with the European Medical Devices Directive 93/42/EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
- Classification:** Class Is.
- Legal manufacturer:** Cedic Srl, Via Liberazione 63/9, IT-20068 Peschiera Borromeo (MI).
- Sustainability:** Environmental management system complies with ISO 14001.
- Materials:** Bag: EVA. Tubing: PVC. All materials are DEHP, BPA and Latex-free*. Additional details on request.
- Country of origin:** Made in Egypt.

*DEHP: Not intentionally added during the manufacturing process. Not formulated with BPA or natural rubber latex.

To obtain further information please contact your Compat® representative or go to www.compat.com
*Unless stated otherwise, all trademarks are used with permission.

AXI00032 REV.2 03 OCTOBER 2025

Compat®